



USER GUIDE

applied
biosystems®
by *life* technologies™

Pathatrix® *Salmonella* spp. Kit (Individual Samples) Linked to PCR or ELISA

For use with the Pathatrix® Auto Instrument

Catalog Number APS50

Publication Number MAN0007089

Revision 11 October 2012

For testing of Food and Environmental samples only.

life
technologies™

The information in this guide is subject to change without notice.

DISCLAIMER

LIFE TECHNOLOGIES CORPORATION AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

Limited Use Label License: Environmental testing, quality control/quality assurance testing, food and agricultural testing

The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product (a) to perform internal research for the sole benefit of the purchaser; and (b) for environmental testing, quality control/quality assurance testing, food and agricultural testing, including reporting results of purchaser's activities in environmental testing, quality control/quality assurance testing, food and agricultural testing for a fee or other commercial consideration. No other right is hereby granted expressly, by implication, or by estoppel. This product is for environmental testing, quality control/ quality assurance testing, food and agricultural testing and research purposes only.

The purchase of this product does not grant the purchaser any additional rights, including (without limitation) the right to transfer or resell the product in any form or the right to use the product as a therapeutic agent or diagnostics test component. For information on obtaining additional rights, please contact outlicensing@lifetech.com or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008.

TRADEMARKS

The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners in the United States and other countries.

BAX is a registered trademark DuPont or its affiliates. R.A.P.I.D. is a registered trademark of Idaho Technology, Inc. Stomacher is a registered trademark of Seward Limited, UK. TECRA and 3M are trademarks of 3M Company Corporation. VIDAS is a registered Trademark of bioMérieux. Whirl-Pak is a registered trademark of Wheatherby/Nasco, Inc.

© 2012 Life Technologies Corporation. All rights reserved.

Contents

Product Information	7
About the kit	7
Kit contents	8
Materials not included in the kit	8
■ CHAPTER 1 Pathatrix® <i>Salmonella</i> spp. Kit Linked to the DuPont® BAX® PCR System	11
Workflow	11
Procedural guidelines	11
Sample enrichment	12
Sample preparation	12
Samples and consumable loading	13
Sample unloading	15
Detection – BAX® PCR protocol	16
Detection – direct plating (optional)	17
Test result interpretation and classification	17
■ CHAPTER 2 Pathatrix® <i>Salmonella</i> spp. Kit Linked to the Idaho Technology R.A.P.I.D.® LT PCR System	19
Workflow	19
Procedural guidelines	19
Sample enrichment	20
Sample preparation	20
Samples and consumable loading	21
Sample unloading	23
Detection – R.A.P.I.D.® LT PCR protocol	24
Detection – direct plating (optional)	25
Test result interpretation and classification	25

■	CHAPTER 3	Pathatrix® <i>Salmonella</i> spp. Kit Linked to the RayAI <i>Salmonella</i> OPTIMA System	27
	Workflow		27
	Procedural guidelines		27
	Sample enrichment		28
	Sample preparation		28
	Samples and consumable loading		29
	Sample unloading		31
	Detection – <i>Salmonella</i> OPTIMA protocol		32
	Detection – direct plating (optional)		33
	Test result interpretation and classification		33
■	CHAPTER 4	Pathatrix® <i>Salmonella</i> spp. Kit Linked to the 3M™ TECRA™ <i>Salmonella</i> Visual Immunoassay	35
	Workflow		35
	Procedural guidelines		35
	Sample enrichment		36
	Sample preparation		36
	Samples and consumable loading		37
	Sample unloading		39
	Detection – TECRA™ <i>Salmonella</i> visual immunoassay protocol		40
	Detection – direct plating (optional)		41
	Test result interpretation and classification		41
■	CHAPTER 5	Pathatrix® <i>Salmonella</i> spp. Kit Linked to the bioMérieux VIDAS® Immunoassay Testing System	43
	Workflow		43
	Procedural guidelines		43
	Sample enrichment		44
	Sample preparation		44
	Samples and consumable loading		45
	Sample unloading		47
	Detection – VIDAS® immunoassay testing protocol		48
	Detection – direct plating (optional)		49
	Test result interpretation and classification		49

■	APPENDIX A	Alternative Enrichment Methods	50
		Specific recommendations for high background microflora	50
		TT broth information	50
		Specific recommendations for milk powder, chocolate, and cocoa-based samples	51
		Specific recommendations for potentially acidic/alkaline samples	51
■	APPENDIX B	Background	52
		Product overview	52
		Description of target microorganisms	52
		Audience	52
		Sampling protocol	52
		Kit sensitivity	52
		Operating conditions	53
■	APPENDIX C	Ordering Information	54
		Related materials from Life Technologies	54
■	APPENDIX D	Safety	55
		Chemical safety	55
		Specific chemical handling	56
		Biological hazard safety	56
		Documentation and Support	57
		Obtaining SDSs	57
		Obtaining Certificates of Analysis	57
		Obtaining support	57
		Limited product warranty	57
		References	58

Product Information

IMPORTANT! Before using this product, read and understand the information in [Appendix D, “Safety” on page 55](#).

CAUTION! *Salmonella* spp. is a Biosafety Level 2 (BSL-2) organism (excluding *S. typhi* and *S. paratyphi*, which are both Biosafety Level 3 [BSL-3]). Care must be taken when handling samples that may contain salmonellae. Laboratory personnel must be adequately trained to handle pathogens before being permitted to analyze samples for *Salmonella* spp. Laboratory personnel must wear appropriate protective equipment, which includes but is not limited to: protective eyewear, face shield, clothing/lab coat, and gloves. Waste should be disposed of in compliance with local and national legislation as appropriate.

About the kit

The Pathatrix® *Salmonella* spp. Kit provides a sample preparation method for presence/absence testing based on the detection of as few as 1–10 cfu (colony forming units)/25–325 g of food sample. An AOAC-validated protocol using this kit for sample preparation followed by detection of microorganisms by selective agar plates can be downloaded from www.lifetechnologies.com/pathatrix (Pathatrix® *Salmonella* spp. Kit (Individual Samples) Linked to Selective Agar Plates, Pub. no. MAN0006974).

This user guide provides 5 alternative protocols for using the Pathatrix® *Salmonella* spp. Kit followed by detection of microorganisms by PCR or ELISA. Presumptive results can be obtained within the following times:

- 20–22 hours, when linked to the DuPont® BAX® PCR system ([Chapter 1](#))
- 17–19 hours, when linked to the Idaho Technology R.A.P.I.D.® LT PCR system ([Chapter 2](#))
- 22–28 hours, when linked to the RayAI *Salmonella* OPTIMA ([Chapter 3](#))
- 22–28 hours, when linked to the 3M™ TECRA™ *Salmonella* visual immunoassay ([Chapter 4](#))
- 22–28 hours, when linked to the bioMérieux VIDAS® immunoassay testing system ([Chapter 5](#))

A presumptive positive isolate should be subsequently confirmed by the use of subculture, as well as appropriate biochemical and serological tests as required.

Once confirmed, the results are reported as:

- *Salmonella* spp. Detected in 25–325 g (sample matrices)
- *Salmonella* spp. Not detected in 25–325 g (sample matrices)

See [Appendix B on page 52](#) for additional background information.

Kit contents

The Pathatrix® *Salmonella* spp. Kit (Cat. no. APS50) contains enough consumable components and Pathatrix® paramagnetic beads to process 50 samples.

Item	Quantity or volume	Storage
Pre-sterilized Sample and Elution Vessel Packs	50 each	Room temperature
Pre-sterilized Capture Phase Packs	50 each	Room temperature
Pre-sterilized Flat Cap Lids	50 each	Room temperature
Anti- <i>Salmonella</i> spp. Antibody-Coated Paramagnetic Beads [†]	2.5 mL (50 tests)	5 ±3°C

[†] The beads have a shelf life of 12 months and are labeled with an expiration date accordingly.

IMPORTANT! Never freeze the Pathatrix® paramagnetic bead suspension. Beads that have been subjected to freezing temperatures may be rendered inactive.

Note: Parts may ship separately depending on configuration and storage conditions.

Materials not included in the kit

The following table includes materials and equipment for using (but not included in) the Pathatrix® consumables kits. Unless otherwise indicated, many of the listed items are available from major laboratory suppliers (MLS).

Item	Source
Equipment	
Incubator, 37 ±1°C	MLS
Pathatrix® Auto Instrument	Life Technologies Cat. no. PATHATRIXAUTO
For use with the BAX® PCR protocol Magnetic Capture Plate	Life Technologies Cat. no. MAGNETICPLATE
For use with the R.A.P.I.D.® LT PCR protocol and VIDAS® protocol DynaMag™-2 Magnet	Life Technologies Cat. no. 123.21D
Forceps, scissors, spatula, knife, and/or scalpel	MLS
Consumables	
Optional for environmental samples: Swabs or sponges	MLS
Sterile bags for enrichment (Whirl-Pak® or Stomacher® bag, or equivalent)	Nasco # B01196WA, Seward product code BA6041, or equivalent

Item	Source
Optional for high-particulate or high-fat-content samples: Sterile filter bags for enrichment (Whirl-Pak® or Stomacher® bag, or equivalent)	Nasco # B01348WA, Seward product code BA6041/STR, or equivalent
Microcentrifuge tubes, PCR clean, 1.5-mL	MLS
Sterile 10-µL disposable loops	MLS
Media	
See Sample Enrichment sections in the protocols (page 12, 20, 28, 36, or 44) or Appendix A on page 50 for recommendations about enrichment media choice. The media is supplied by several manufacturers (e.g., Oxoid [product codes shown], Difco, and Merck) in a dehydrated form and should be prepared according to the manufacturer's instructions.	
Buffered peptone water	Oxoid product code CM0509
For samples with high background microflora: Tetrathionate (TT) broth base	Oxoid product code CM0029
For milk powder, chocolate, or cocoa-based samples: UHT skim milk Brilliant green (CAS 633-03-4)	Food retail store MLS
Selective agar	
XLD agar	Oxoid product code CM0469
Brilliant green agar (modified)	Oxoid product code CM0329
Reagents	
PBS, 10X, pH 7.4	Life Technologies Cat. no. AM9624 or AM9625


Product Information

Materials not included in the kit

1

Pathatrix[®] *Salmonella* spp. Kit Linked to the DuPont[®] BAX[®] PCR System

Workflow

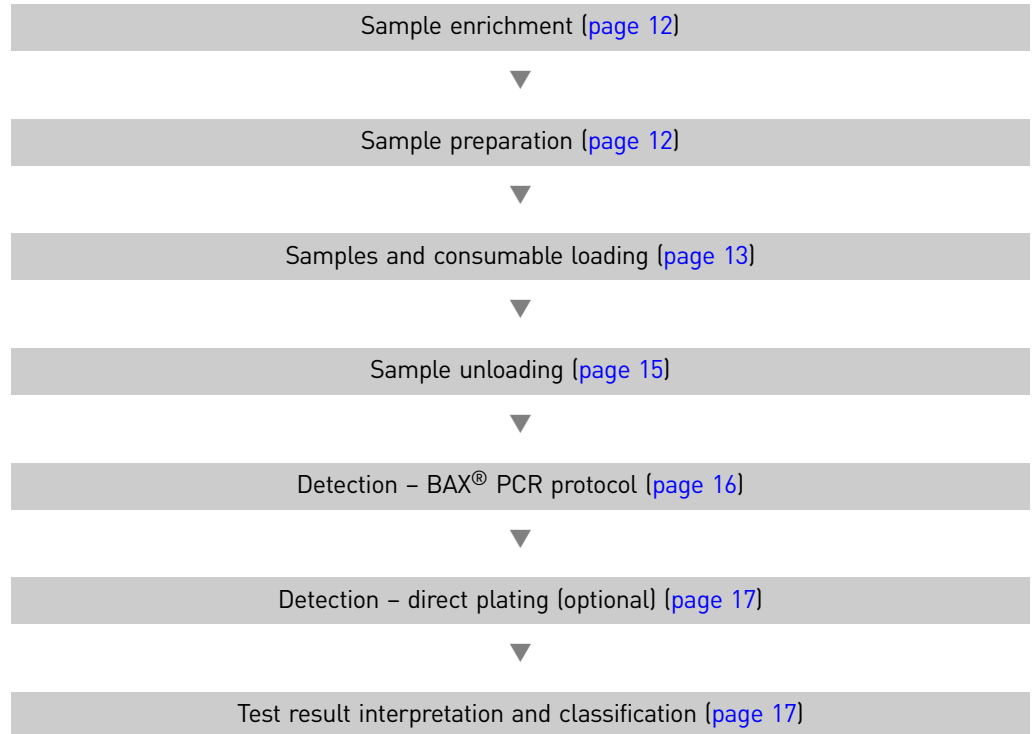

Procedural guidelines

- Use aseptic technique and good laboratory practices at all times.
- A facemask should be worn when weighing out powders.
- Care must be taken when boiling agar prior to autoclaving, and heat-resistant gloves should be worn when handling hot flasks of liquid.
- Take care when handling plates or vessels that may contain microorganisms.
- Avoid generating aerosols, as pathogenic organisms may be present.
- Used or unused reagents, used media, sample enrichments, as well as other contaminated disposable materials should be disposed of following procedures for infectious or potentially infectious products.
- Sample enrichments might be contaminated with pathogenic organisms infectious to humans, so all waste must be treated as biohazardous and handled and disposed using safe laboratory practices, in accordance and compliance with all appropriate regulations.

Sample enrichment

Note: Certain food types and swabs/sponges can benefit from an alternative enrichment strategy (see [Appendix A, “Alternative Enrichment Methods”](#)).

1. Prepare a 1:10 dilution of the food sample in the appropriate **prewarmed** ($37 \pm 1^\circ\text{C}$) enrichment media in a sterile bag.

Note: For example, add 25 g of food sample to 225 mL of prewarmed media or add 325 g of food sample to 2925 mL of prewarmed media.

IMPORTANT! It is critical that the enrichment media is **prewarmed** to $37 \pm 1^\circ\text{C}$ prior to adding the food sample. To prevent cooling, all samples should then immediately be placed in a 37°C incubator.

2. Homogenize the sample by massaging the sample between the fingers for 10–15 seconds to disperse any large clumps of material (Narang and Cray, 2006).
3. Incubate at $37 \pm 1^\circ\text{C}$ for a **minimum** of **16 hours**.

Sample preparation

1. Remove the Sample Vessel and Elution Vessel from the consumable kit packaging and place into the Sample Vessel Holder.
2. Partially remove the lids from both vessels, making a large enough opening to allow the sample and wash buffer to be dispensed into the vessels.
3. Place 50 μL of your sample in the Sample Vessel.


Note: If the samples are highly particulate or has a high fat content, a Seward plain sterile bag with internal strainer may be used (Seward Product Code BA6041/STR).

4. Store the individual enriched samples at $5 \pm 3^\circ\text{C}$ for potential reanalysis until the test result is confirmed.


Note: Do not store for more than 32 hours. If refrigerated after pre-enrichment, samples should be rewarmed to $37 \pm 1^\circ\text{C}$ prior to removal of aliquots for analysis.

Samples and consumable loading


1. Add PBS (Cat. no. AM9624, diluted to 1X) to the fill line of the Elution Vessel (approximately 35 mL of 1X PBS).
2. Place the lids on the Sample and Elution Vessels, making sure that the vessels are sealed all the way around.
3. Ensure the Pathatrix® paramagnetic beads (included in the kit) are fully resuspended by agitating the bead vial (for example, by vortexing bead vial or inversion of sealed vial), and add 50 µL of the Pathatrix® paramagnetic bead suspension into the spout on the lid of the Sample Vessel.

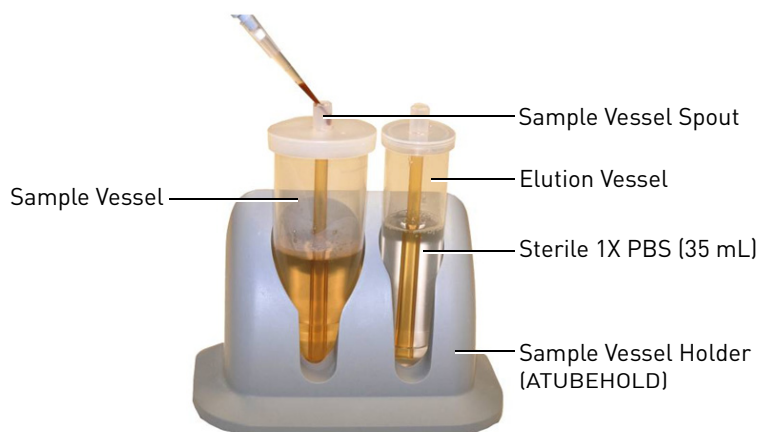

4. Remove the capture-phase kit from the bag and orient with the valve plunger pointing left. Connect the valve firmly to the lids of the Sample and Elution Vessels.

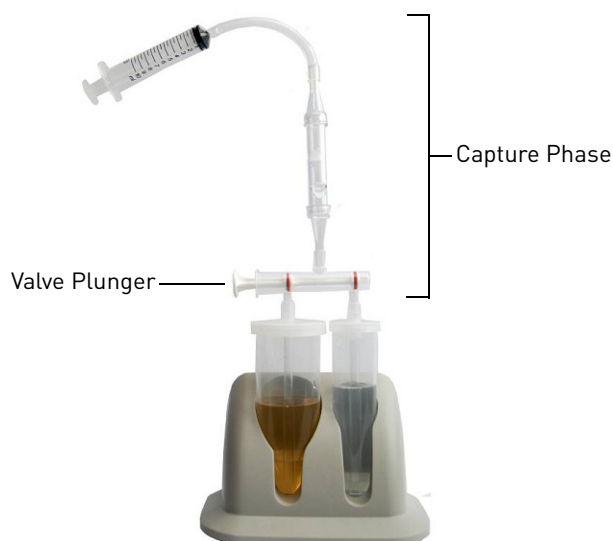

5. Holding both vessels, lift the assembled vessels and attached capture-phase kit out of the Sample Vessel Holder.
6. Place the vessels into the Cartridge, pushing them firmly in place from the bottom upwards.
7. Firmly push the remainder of the kit into the Cartridge, ensuring that the valve, the capture phase, and the syringe are all held securely in the molded recess of the Cartridge.
8. Push the magnet slider across into the locking position and press the magnet release button on the side of the Cartridge to ensure the kit is positioned correctly and the magnets can freely disengage from the capture phase.

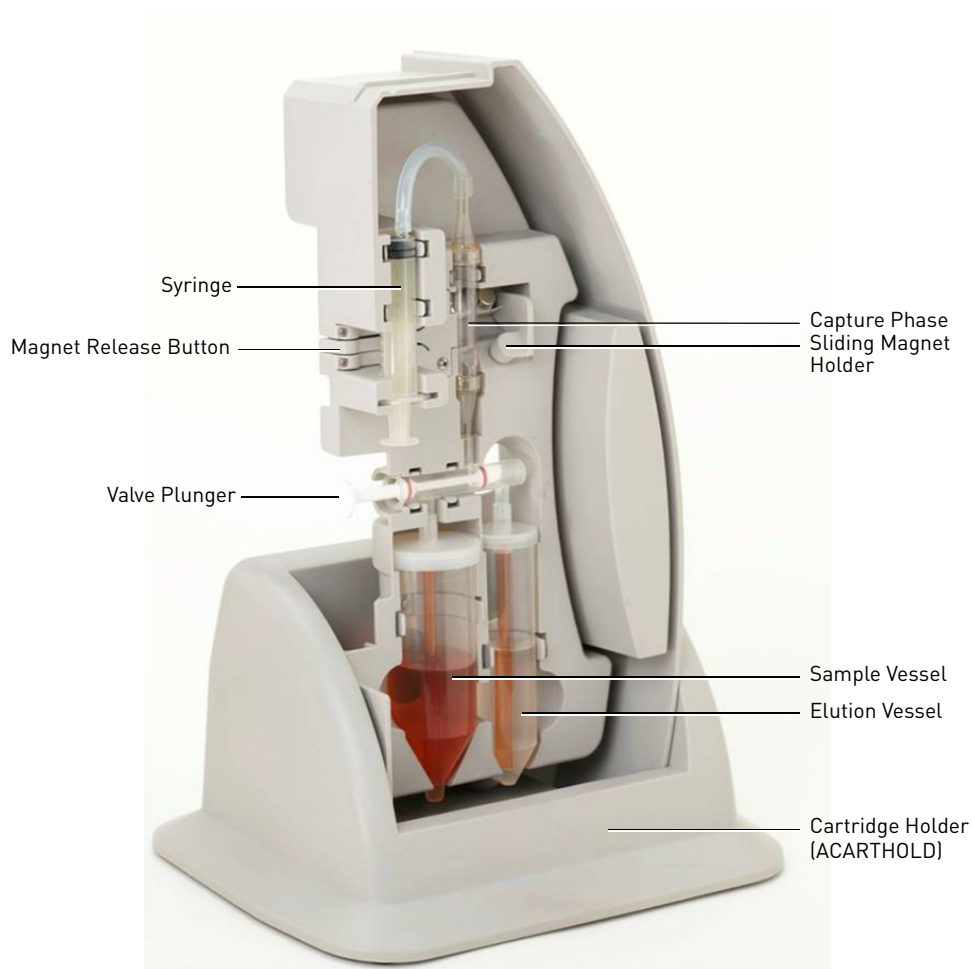

9. Reset the magnets into the locking position.
10. Insert the Cartridge into the Pathatrix® Auto Instrument until it clicks into the locking position.
11. Press the numbered button above the appropriate Cartridge to start the run (approximately 14 minutes). The associated LED will turn green to indicate the run has started.

Sample unloading


1. At the end of the run, the LED will flash red and green alternately.
2. Press the button above the appropriate Cartridge to initialize the draining step (approximately 1 minute).
3. When the draining step is complete and the LED is illuminated red, remove the Cartridge by pulling it out, away from the instrument.
4. Remove the syringe from the Cartridge and carefully pull out the rest of the kit, starting at the top and working downwards. When removing the Sample and Elution Vessels from the Cartridge, hold firmly by the vessels themselves to prevent spillage.
5. Place both vessels in the Tube Rack (also known as the Sample Vessel Holder or rack).
6. Remove the lid from the Elution Vessel and, leaving the Elution Vessel in place in the rack, lift away the rest of the consumable, including the Sample Vessel, and discard.
7. Place the Elution Vessel in the rack and cap with a flat lid (provided in the kit). Leave in the rack for 1 minute to allow capture of the Pathatrix® paramagnetic beads.
8. Remove all the liquid from the Elution Vessel, without removing the vessel from the rack, taking care not to disturb the captured Pathatrix® paramagnetic beads.
9. Remove the Elution Vessel from the vessel holder, add 100 µL of PBS into the Elution Vessel, and resuspend the Pathatrix® paramagnetic beads.
10. Appropriate aliquots of the Pathatrix® paramagnetic bead suspension can then be immediately analyzed using the laboratory's chosen pathogen detection method.

Note: The Pathatrix® paramagnetic bead suspension may be retained for later testing if necessary. The beads should be stored in sterile 1.5-mL microcentrifuge tubes AWAY from magnets (for example, away from Pathatrix® vessel holders) at 5 ±3°C for up to 24 hours.

Detection – BAX® PCR protocol

IMPORTANT! It is solely the operator's responsibility to refer to the package insert instructions and instrument user manual for information regarding the correct set up and operation of the BAX® PCR system.

1. Pipet 25 µL of the Pathatrix® paramagnetic bead suspension into BAX® lysis tubes and proceed according to the standard BAX® lysis procedure.
2. Once the lysis step is complete, immediately place the rack and tubes onto the Magnetic Capture Plate (Cat. no. MAGNETICPLATE) and leave for at least 1 minute to allow the Pathatrix® paramagnetic beads to accumulate on the bottom of the tube.
3. Transfer only **bead-free lysate** into the BAX® PCR reaction tubes.
Note: Target DNA, if present, will be in the **bead-free** supernatant.
4. To proceed, refer to the DuPont® Qualicon operating instructions.

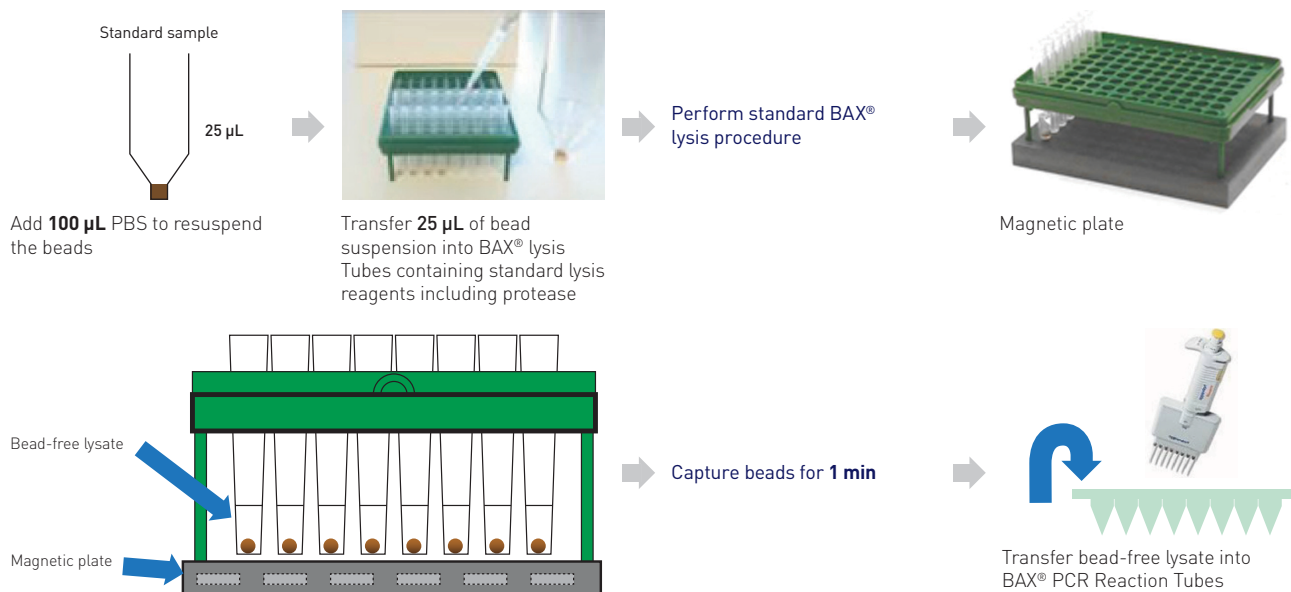

If a positive PCR result is obtained, an aliquot of the Pathatrix® paramagnetic beads should be plated out (see the following section, “Detection – direct plating (optional)”).

Detection – direct plating (optional)

Note: We recommend streaking the sample out to generate individual colonies, as opposed to spread plating.

1. Using an inoculation loop, streak 10 µL of the unlysed Pathatrix® paramagnetic bead suspension over a well-dried XLD (xylose lysine desoxycholate) agar plate and another 10 µL onto an appropriate second selective plate medium.

Note: The laboratory may choose which medium to use, but the second selective plate should be any other solid selective medium complimentary to XLD and especially appropriate for the isolation of lactose-positive *Salmonella*, *Salmonella typhi*, and *Salmonella paratyphi* strains.

2. Allow the plates to dry for approximately 10 minutes then invert and incubate at the required temperature for 18–24 hours or as recommended by the manufacturer.

A presumptive positive result is defined as the isolation of typical, suspicious, or atypical *Salmonella* spp. colonies on the agar plates used. A presumptive positive isolate should be subsequently confirmed by the use of subculture and appropriate biochemical and serological tests as detailed in USDA Microbiology Laboratory Guidebook (MLG) 4.04 for cooked ham or FDA Bacteriological Analytical Manual (BAM), Chapter 5, sections D and E for tomatoes and chocolate (see [“References” on page 58](#)).

Test result interpretation and classification

The Pathatrix® *Salmonella* spp. Kit is designed as a sample preparation method for presence/absence detection of *Salmonella* spp. in food matrices

Using the Pathatrix® *Salmonella* spp. Kit linked to the BAX® PCR system, presumptive results can be obtained, prior to confirmation, within 20–22 hours.


Once confirmed, the results are reported as:

- *Salmonella* spp. **Detected** in 25–325 g (sample matrices)
- *Salmonella* spp. **Not detected** in 25–325 g (sample matrices)

2

Pathatrix® *Salmonella* spp. Kit Linked to the Idaho Technology R.A.P.I.D.® LT PCR System

Workflow

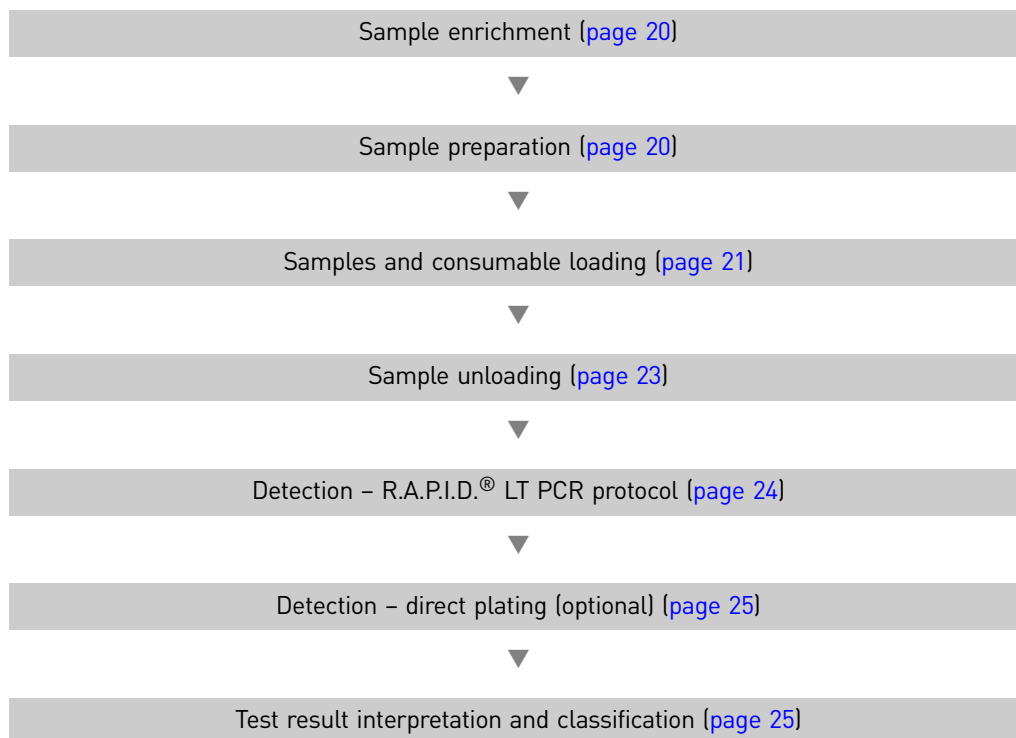

Procedural guidelines

- Use aseptic technique and good laboratory practices at all times.
- A facemask should be worn when weighing out powders.
- Care must be taken when boiling agar prior to autoclaving, and heat-resistant gloves should be worn when handling hot flasks of liquid.
- Take care when handling plates or vessels that may contain microorganisms.
- Avoid generating aerosols, as pathogenic organisms may be present.
- Used or unused reagents, used media, sample enrichments, as well as other contaminated disposable materials should be disposed of following procedures for infectious or potentially infectious products.
- Sample enrichments might be contaminated with pathogenic organisms infectious to humans, so all waste must be treated as biohazardous and handled and disposed using safe laboratory practices, in accordance and compliance with all appropriate regulations.

Sample enrichment

Note: Certain food types and swabs/sponges can benefit from an alternative enrichment strategy (see [Appendix A, “Alternative Enrichment Methods”](#)).

1. Prepare a 1:10 dilution of the food sample in the appropriate **prewarmed** ($37 \pm 1^\circ\text{C}$) enrichment media in a Stomacher® bag.

Note: For example, add 25 g of food sample to 225 mL of prewarmed media or add 325 g of food sample to 2925 mL of prewarmed media.

IMPORTANT! It is critical that the enrichment media is **prewarmed** to $37 \pm 1^\circ\text{C}$ prior to adding the food sample. To prevent cooling, all samples should then immediately be placed in a 37°C incubator.

2. Homogenize the sample by massaging the sample between the fingers for 10–15 seconds to disperse any large clumps of material (Narang and Cray, 2006).
3. Incubate at $37 \pm 1^\circ\text{C}$ for a **minimum** of **16 hours**.

Sample preparation

1. Remove the Sample Vessel and Elution Vessel from the consumable kit packaging and place into the Sample Vessel Holder.
2. Partially remove the lids from both vessels, making a large enough opening to allow the sample and wash buffer to be dispensed into the vessels.
3. Place 50 mL of your sample in the Sample Vessel.


Note: If the samples are highly particulate or has a high fat content, a Seward plain sterile bag with internal strainer may be used (Seward Product Code BA6041/STR).

4. Store the individual enriched samples at $5 \pm 3^\circ\text{C}$ for potential reanalysis until the test result is confirmed.


Note: Do not store for more than 32 hours. If refrigerated after pre-enrichment, samples should be re-warmed to $37 \pm 1^\circ\text{C}$ prior to removal of aliquots for analysis.

Samples and consumable loading


1. Add PBS (Cat. no. AM9624, diluted to 1X) to the fill line of the Elution Vessel (approximately 35 mL of 1X PBS).
2. Place the lids on the Sample and Elution Vessels, making sure that the vessels are sealed all the way around.
3. Ensure the Pathatrix® paramagnetic beads (included in the kit) are fully resuspended by agitating the bead vial (for example, by vortexing bead vial or inversion of sealed vial), and add 50 µL of the Pathatrix® paramagnetic bead suspension into the spout on the lid of the Sample Vessel.

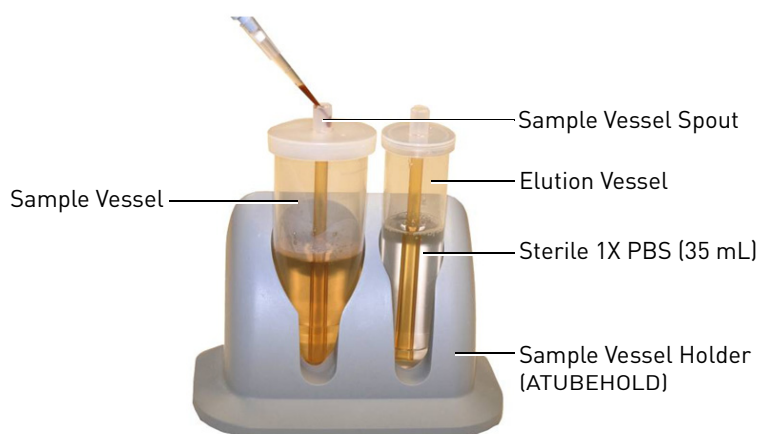

4. Remove the capture-phase kit from the bag and orient with the valve plunger pointing left. Connect the valve firmly to the lids of the Sample and Elution Vessels.

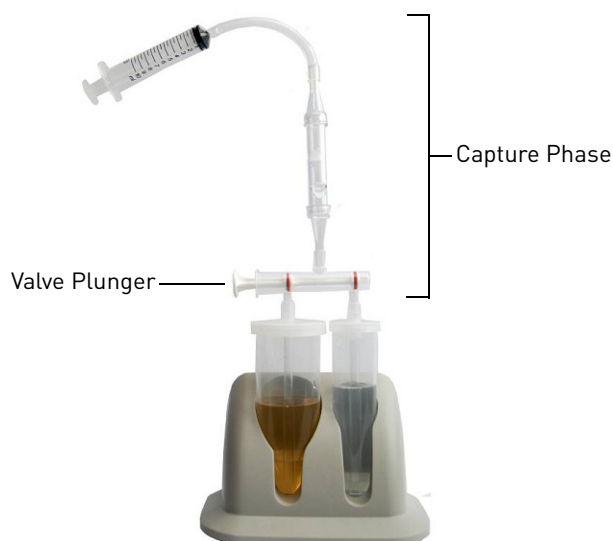

5. Holding both vessels, lift the assembled vessels and attached capture-phase kit out of the Sample Vessel Holder.
6. Place the vessels into the Cartridge, pushing them firmly in place from the bottom upwards.
7. Firmly push the remainder of the kit into the Cartridge, ensuring that the valve, the capture phase, and the syringe are all held securely in the molded recess of the Cartridge.
8. Push the magnet slider across into the locking position and press the magnet release button on the side of the Cartridge to ensure the kit is positioned correctly and the magnets can freely disengage from the capture phase.

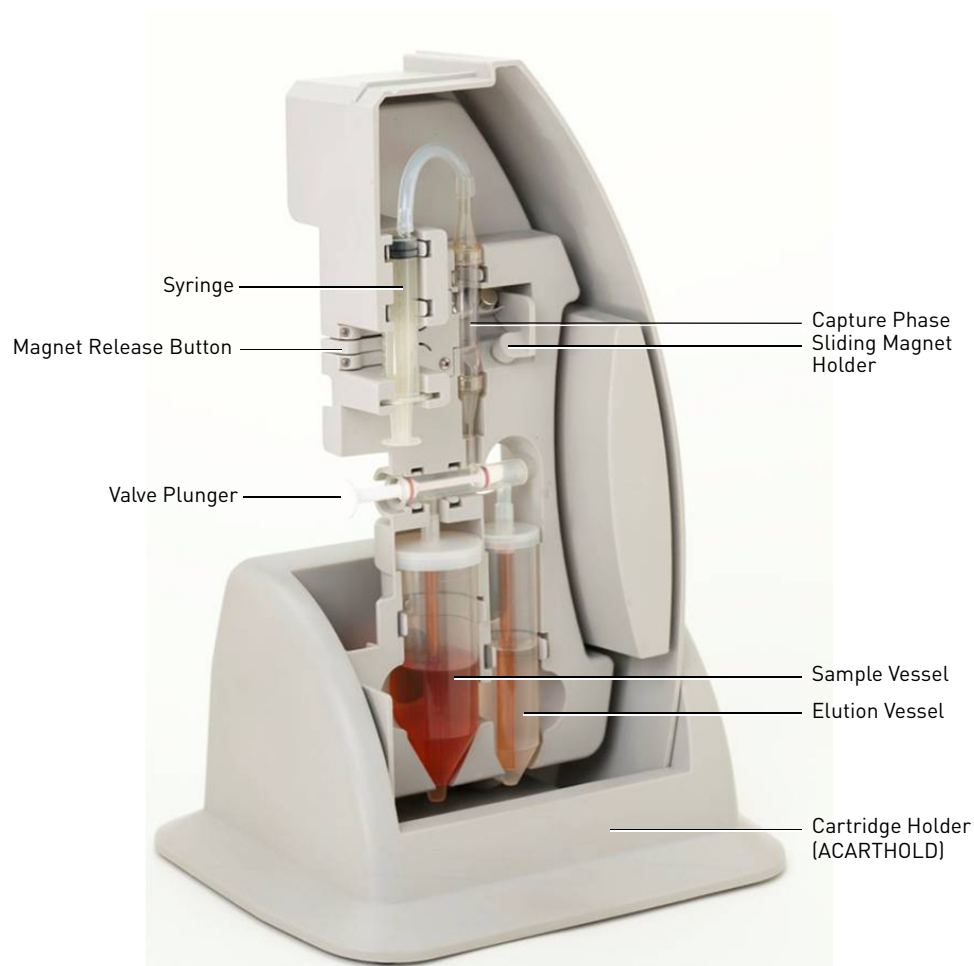

9. Reset the magnets into the locking position.
10. Insert the Cartridge into the Pathatrix® Auto Instrument until it clicks into the locking position.
11. Press the numbered button above the appropriate Cartridge to start the run (approximately 14 minutes). The associated LED will turn green to indicate the run has started.

Sample unloading


1. At the end of the run, the LED will flash red and green alternately.
2. Press the button above the appropriate Cartridge to initialize the draining step (approximately 1 minute).
3. When the draining step is complete and the LED is illuminated red, remove the Cartridge by pulling it out, away from the instrument.
4. Remove the syringe from the Cartridge and carefully pull out the rest of the kit, starting at the top and working downwards. When removing the Sample and Elution Vessels from the Cartridge, hold firmly by the vessels themselves to prevent spillage.
5. Place both vessels in the Tube Rack (also known as the Sample Vessel Holder or rack).
6. Remove the lid from the Elution Vessel and, leaving the Elution Vessel in place in the rack, lift away the rest of the consumable, including the Sample Vessel, and discard.
7. Place the Elution Vessel in the rack and cap with a flat lid (provided in the kit). Leave in the rack for 1 minute to allow capture of the Pathatrix® paramagnetic beads.
8. Remove all the liquid from the Elution Vessel, without removing the vessel from the rack, taking care not to disturb the captured Pathatrix® paramagnetic beads.
9. Remove the Elution Vessel from the vessel holder, add 100 µL of PBS into the Elution Vessel, and resuspend the Pathatrix® paramagnetic beads.
10. Appropriate aliquots of the Pathatrix® paramagnetic bead suspension can then be immediately analyzed using the laboratory's chosen pathogen detection method.

Note: The Pathatrix® paramagnetic bead suspension may be retained for later testing if necessary. The beads should be stored in sterile 1.5-mL microcentrifuge tubes AWAY from magnets (for example, away from Pathatrix® vessel holders) at 5 ±3°C for up to 24 hours.

Detection – R.A.P.I.D.® LT PCR protocol

IMPORTANT! It is solely the operator's responsibility to refer to the package insert instructions and instrument user manual for information regarding the correct set up and operation of the R.A.P.I.D.® LT PCR system.

1. Add 25 μ L of the resuspended Pathatrix® paramagnetic beads into a "Blue Capped Bead Tube" and vortex on the Disruptor Genie® vortexer (equipped with the TurboMix™ attachment) for **5 minutes** on the highest setting.
2. Once the bead-beating lysis step is complete, immediately place the Bead Tubes into the DynaMag™-2 Magnet (Cat. no. 123.21D).
3. Wait for at least 1 minute to allow the Pathatrix® paramagnetic bead debris to be drawn out of suspension, thereby producing the **bead-free** supernatant.
Note: If present, target DNA will be in the bead-free lysate.
4. Pipet 10 μ L of Reconstitution Buffer into the "Unknown" test vial.
Note: Once pellets have been reconstituted, samples should be processed immediately.
5. Pipet 10 μ L of bead-free supernatant into the PCR "Unknown" test reaction vial. Mix well to fully resuspend the PCR reagent pellet.
Note: Take care not to transfer any glass beads from the Blue Capped Bead-Beating Tubes into the PCR reagent vial.
6. Transfer 18 μ L of the reagent/sample suspension, and add into the top of a glass PCR capillary tube. Using the "Capping Tool" provided, place the cap on the sample capillary tube and push down firmly.
7. To proceed, refer to the R.A.P.I.D.® LT PCR operating instructions.

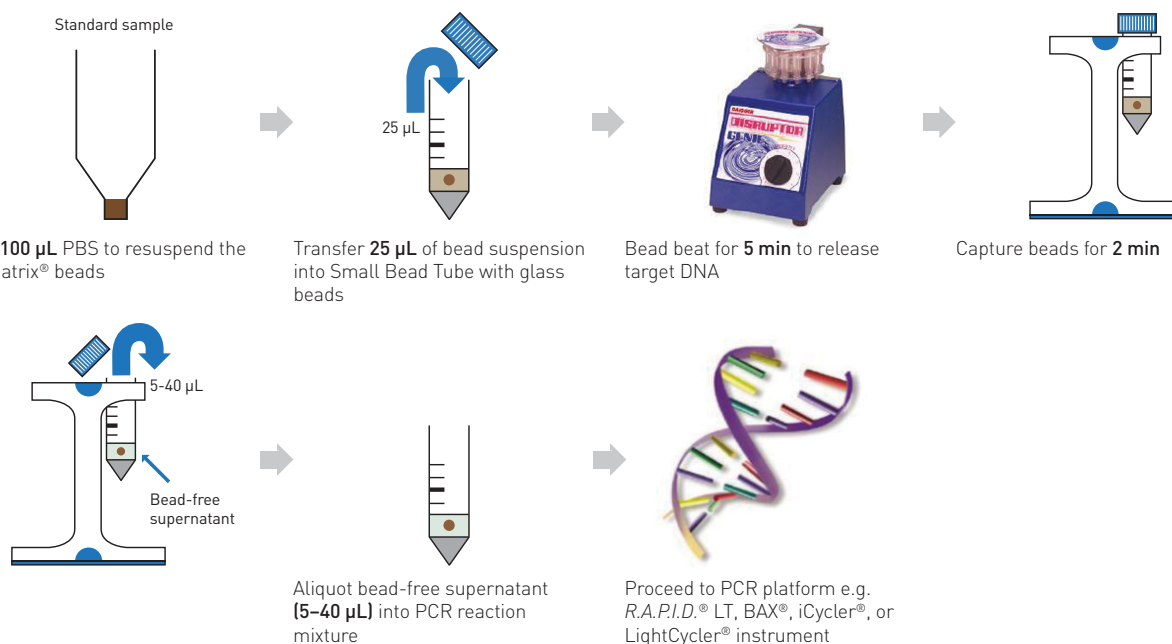

Review the amplification curves and melt peaks as described in the R.A.P.I.D.® LT PCR software “Screen Shot Protocol.” If any Amplification Curve exhibits an increase in fluorescence and/or displays a potential Melt Peak in the absence of a positive PCR software determination, a repeat PCR analysis should be carried out by repeating the R.A.P.I.D.® LT PCR steps in this section (see [page 24](#)).

If a positive PCR result is obtained, an aliquot of the Pathatrix® paramagnetic beads should be plated out (see the following section, “[Detection – direct plating \(optional\)](#)”).

Detection – direct plating (optional)

Note: We recommend streaking the sample out to generate individual colonies, as opposed to spread plating.

1. Using an inoculation loop, streak 10 µL of the unlysed Pathatrix® paramagnetic bead suspension over a well-dried XLD (xylose lysine desoxycholate) agar plate and another 10 µL onto an appropriate second selective plate medium.

Note: The laboratory may choose which medium to use, but the second selective plate should be any other solid selective medium complimentary to XLD and especially appropriate for the isolation of lactose-positive *Salmonella*, *Salmonella typhi*, and *Salmonella paratyphi* strains.

2. Allow the plates to dry for approximately 10 minutes then invert and incubate at the required temperature for 18–24 hours or as recommended by the manufacturer.

A presumptive positive result is defined as the isolation of typical, suspicious, or atypical *Salmonella* spp. colonies on the agar plates used. A presumptive positive isolate should be subsequently confirmed by the use of subculture and appropriate biochemical and serological tests as detailed in USDA Microbiology Laboratory Guidebook (MLG) 4.04 for cooked ham or FDA Bacteriological Analytical Manual (BAM), Chapter 5, sections D and E for tomatoes and chocolate (see “[References](#)” on [page 58](#)).

Test result interpretation and classification

The Pathatrix® *Salmonella* spp. Kit is designed as a sample preparation method for presence/absence detection of *Salmonella* spp. in food matrices

Using the Pathatrix® *Salmonella* spp. Kit linked to the R.A.P.I.D.® LT PCR system, presumptive results can be obtained, prior to confirmation, within 17–19 hours.


Once confirmed, the results are reported as:

- *Salmonella* spp. **Detected** in 25–325 g (sample matrices)
- *Salmonella* spp. **Not detected** in 25–325 g (sample matrices)

3

Pathatrix® *Salmonella* spp. Kit Linked to the RayAl *Salmonella* OPTIMA System

Workflow

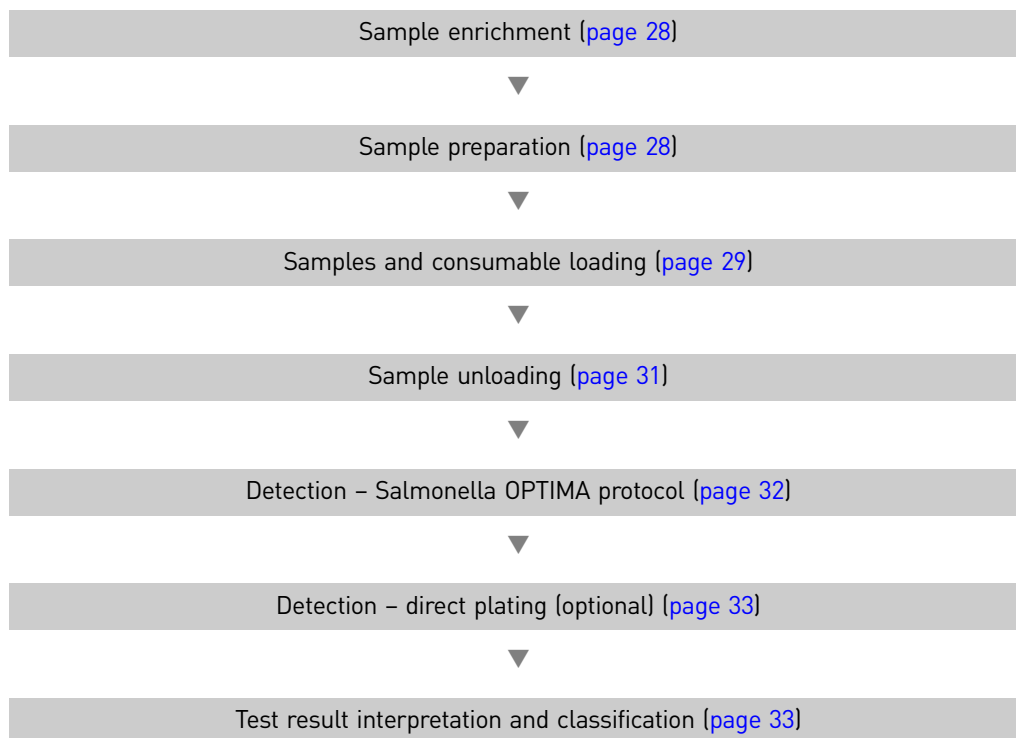

Procedural guidelines

- Use aseptic technique and good laboratory practices at all times.
- A facemask should be worn when weighing out powders.
- Care must be taken when boiling agar prior to autoclaving, and heat-resistant gloves should be worn when handling hot flasks of liquid.
- Take care when handling plates or vessels that may contain microorganisms.
- Avoid generating aerosols, as pathogenic organisms may be present.
- Used or unused reagents, used media, sample enrichments, as well as other contaminated disposable materials should be disposed of following procedures for infectious or potentially infectious products.
- Sample enrichments might be contaminated with pathogenic organisms infectious to humans, so all waste must be treated as biohazardous and handled and disposed using safe laboratory practices, in accordance and compliance with all appropriate regulations.

Sample enrichment

Note: Certain food types and swabs/sponges can benefit from an alternative enrichment strategy (see [Appendix A, “Alternative Enrichment Methods”](#)).

1. Prepare a 1:10 dilution of the food sample in the appropriate **prewarmed** ($37 \pm 1^\circ\text{C}$) enrichment media in a sterile bag.

Note: For example, add 25 g of food sample to 225 mL of prewarmed media or add 325 g of food sample to 2925 mL of prewarmed media.

IMPORTANT! It is critical that the enrichment media is **prewarmed** to $37 \pm 1^\circ\text{C}$ prior to adding the food sample. To prevent cooling, all samples should then immediately be placed in a 37°C incubator.

2. Homogenize the sample by massaging the sample between the fingers for 10–15 seconds to disperse any large clumps of material (Narang and Cray, 2006).
3. Incubate at $37 \pm 1^\circ\text{C}$ for a **minimum** of **18–24 hours**.

Sample preparation

1. Remove the Sample Vessel and Elution Vessel from the consumable kit packaging and place into the Sample Vessel Holder.
2. Partially remove the lids from both vessels, making a large enough opening to allow the sample and wash buffer to be dispensed into the vessels.
3. Place 50 mL of your sample in the Sample Vessel.


Note: If the samples are highly particulate or has a high fat content, a Seward plain sterile bag with internal strainer may be used (Seward Product Code BA6041/STR).

4. Store the individual enriched samples at $5 \pm 3^\circ\text{C}$ for potential reanalysis until the test result is confirmed.


Note: Do not store for more than 32 hours. If refrigerated after pre-enrichment, samples should be rewarmed to $37 \pm 1^\circ\text{C}$ prior to removal of aliquots for analysis.

Samples and consumable loading


1. Add PBS (Cat. no. AM9624, diluted to 1X) to the fill line of the Elution Vessel (approximately 35 mL of 1X PBS).
2. Place the lids on the Sample and Elution Vessels, making sure that the vessels are sealed all the way around.
3. Ensure the Pathatrix® paramagnetic beads (included in the kit) are fully resuspended by agitating the bead vial (for example, by vortexing bead vial or inversion of sealed vial), and add 50 µL of the Pathatrix® paramagnetic bead suspension into the spout on the lid of the Sample Vessel.

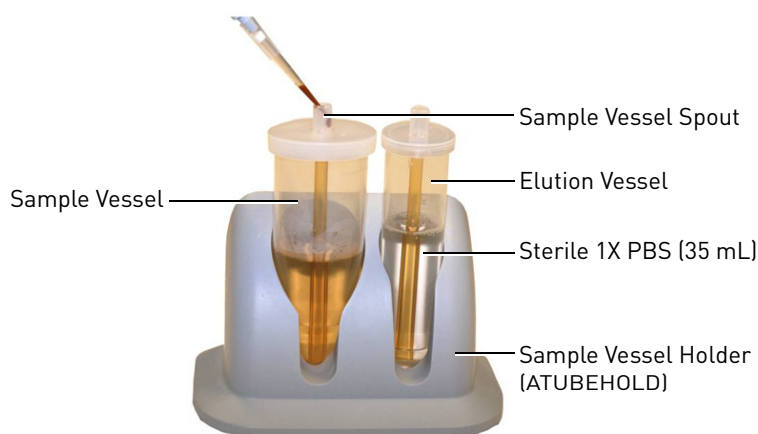

4. Remove the capture-phase kit from the bag and orient with the valve plunger pointing left. Connect the valve firmly to the lids of the Sample and Elution Vessels.

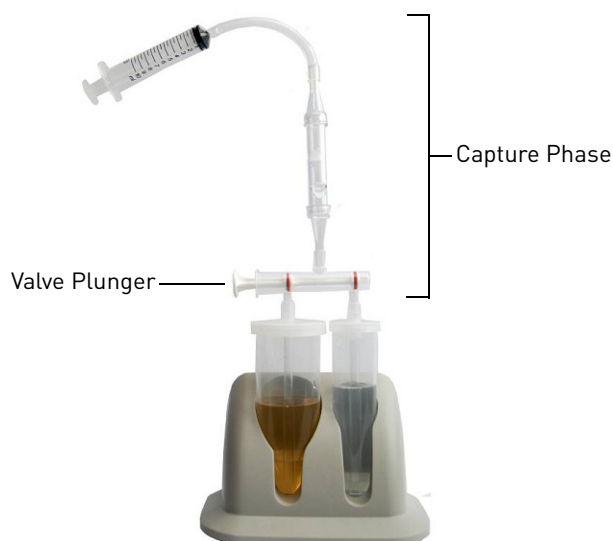

5. Holding both vessels, lift the assembled vessels and attached capture-phase kit out of the Sample Vessel Holder.
6. Place the vessels into the Cartridge, pushing them in firmly from the bottom upwards.
7. Firmly push the remainder of the kit into the Cartridge, ensuring that the valve, the capture phase, and the syringe are all held securely in the molded recess of the Cartridge.
8. Push the magnet slider across into the locking position and press the magnet release button on the side of the Cartridge to ensure the kit is positioned correctly and the magnets can freely disengage from the capture phase.

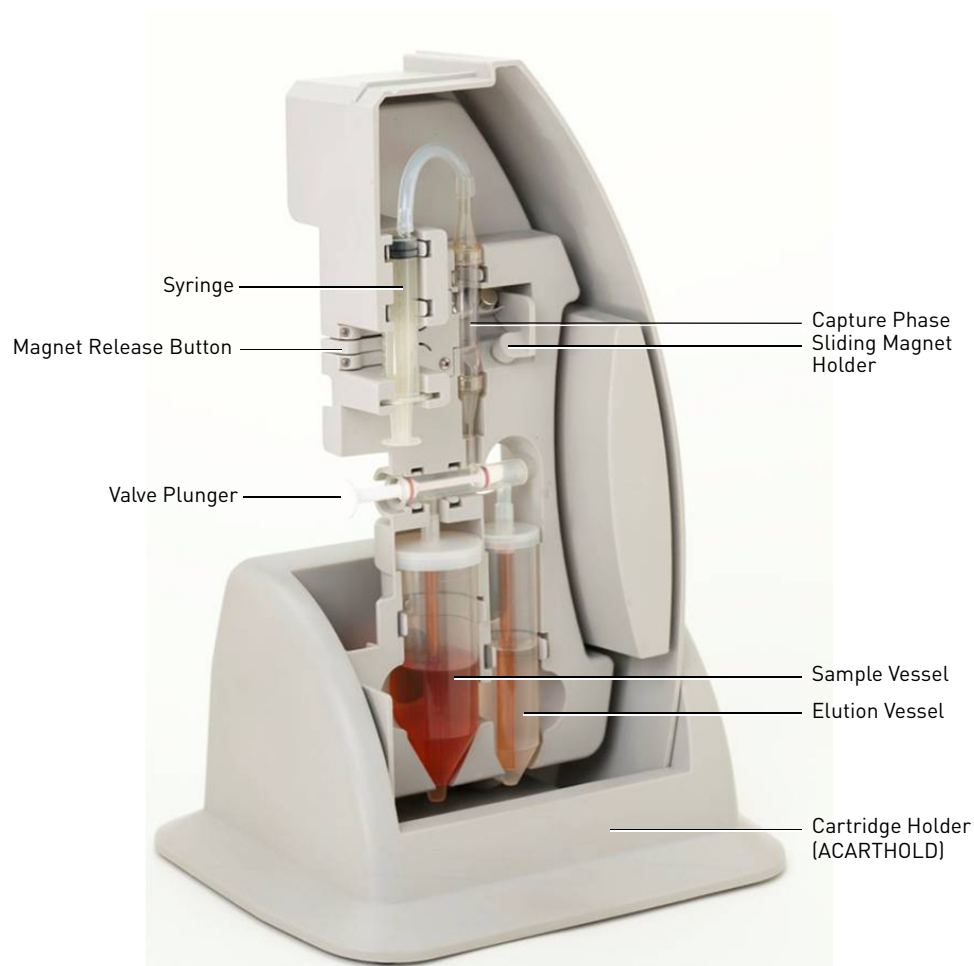

9. Reset the magnets into the locking position.
10. Insert the Cartridge into the Pathatrix® Auto Instrument until it clicks into the locking position.
11. Press the numbered button above the appropriate Cartridge to start the run (approximately 14 minutes). The associated LED will turn green to indicate the run has started.


Sample unloading

1. At the end of the run, the LED will flash red and green alternately.
2. Press the button above the appropriate Cartridge to initialize the draining step (approximately 1 minute).
3. When the draining step is complete and the LED is illuminated red, remove the Cartridge by pulling it out, away from the instrument.
4. Remove the syringe from the Cartridge and carefully pull out the rest of the kit, starting at the top and working downwards. When removing the Sample and Elution Vessels from the Cartridge, hold firmly by the vessels themselves to prevent spillage.
5. Place both vessels in the Tube Rack (also known as the Sample Vessel Holder or rack).
6. Remove the lid from the Elution Vessel and, leaving the Elution Vessel in place in the rack, lift away the rest of the consumable, including the Sample Vessel, and discard.
7. Place the Elution Vessel in the rack and cap with a flat lid (provided in the kit). Leave in the rack for 1 minute to allow capture of the Pathatrix® paramagnetic beads.
8. Remove all the liquid from the Elution Vessel, without removing the vessel from the rack, taking care not to disturb the captured Pathatrix® paramagnetic beads.
9. Remove the Elution Vessel from the vessel holder, use 1200 µL of **pre-warmed** BPW to resuspend the Pathatrix® paramagnetic beads, and transfer 1100 µL to a microcentrifuge tube.
Note: 100 µL of the Pathatrix® paramagnetic beads can be retained for direct plating prior to the additional 2-hour growth step.
10. Vortex the tube for 10 seconds.
11. Incubate the tube for 2 hours at 37 ±1°C.
12. Remove the tube from the incubator and vortex for 10 seconds.

Detection – *Salmonella* OPTIMA protocol

IMPORTANT! It is solely the operator's responsibility to refer to the package insert instructions and instrument user manual for information regarding the correct set up and operation of the *Salmonella* OPTIMA system.

1. Ensure that the tubes are securely capped and then heat in a water bath for 15 minutes at 95°C to release target antigens.
2. After heating, immediately place the tube into the DynaMag™-2 Magnet (Cat. no. 123.21D) for 2 minutes to allow the Pathatrix® paramagnetic beads to be pulled out of the lysate.
3. Once the Pathatrix® paramagnetic beads have aggregated against the magnet, open the tube and transfer 1100 µL of the bead-free supernatant into the designated well of the RayAl microplate.
4. To proceed, refer to the *Salmonella* OPTIMA operating instructions.

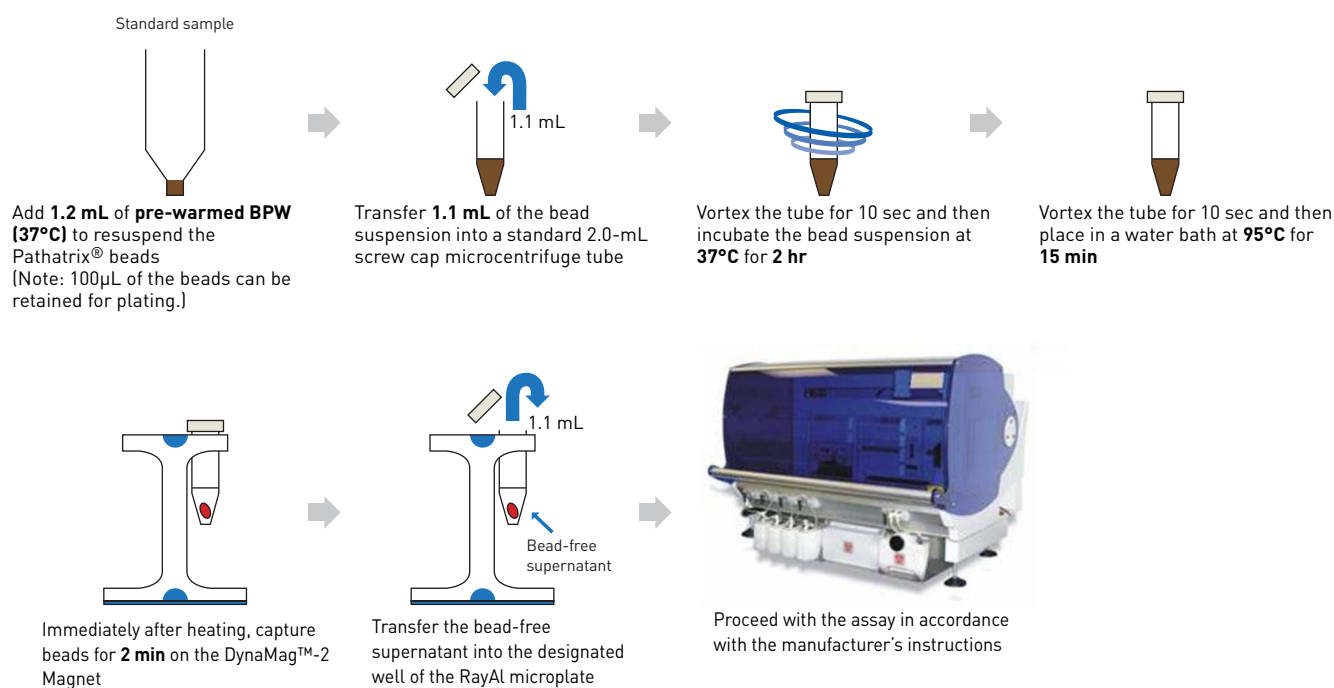

If a positive ELISA result is obtained, an aliquot of the Pathatrix® paramagnetic beads should be plated out (see the following section, “[Detection – direct plating \(optional\)](#)”).

Detection – direct plating (optional)

Note: We recommend streaking the sample out to generate individual colonies, as opposed to spread plating.

1. Streak 50 µL of the unlysed Pathatrix® paramagnetic bead suspension over a well-dried XLD (xylose lysine desoxycholate) agar plate and another 50 µL onto an appropriate second selective plate medium.

Note: The laboratory may choose which medium to use, but the second selective plate should be any other solid selective medium complimentary to XLD and especially appropriate for the isolation of lactose-positive *Salmonella*, *Salmonella typhi*, and *Salmonella paratyphi* strains.

2. Allow the plates to dry for approximately 10 minutes then invert and incubate at the required temperature for 18–24 hours or as recommended by the manufacturer.

A presumptive positive result is defined as the isolation of typical, suspicious, or atypical *Salmonella* spp. colonies on the agar plates used. A presumptive positive isolate should be subsequently confirmed by the use of subculture and appropriate biochemical and serological tests as detailed in USDA Microbiology Laboratory Guidebook (MLG) 4.04 for cooked ham or FDA Bacteriological Analytical Manual (BAM), Chapter 5, sections D and E for tomatoes and chocolate (see [“References” on page 58](#)).

Test result interpretation and classification

The Pathatrix® *Salmonella* spp. Kit is designed as a sample preparation method for presence/absence detection of *Salmonella* spp. in food matrices

Using the Pathatrix® *Salmonella* spp. Kit linked to RayAl *Salmonella* OPTIMA, presumptive results can be obtained, prior to confirmation, within 22–28 hours.


Once confirmed, the results are reported as:

- *Salmonella* spp. **Detected** in 25–325 g (sample matrices)
- *Salmonella* spp. **Not detected** in 25–325 g (sample matrices)

4

Pathatrix® *Salmonella* spp. Kit Linked to the 3M™ TECRA™ *Salmonella* Visual Immunoassay

Workflow

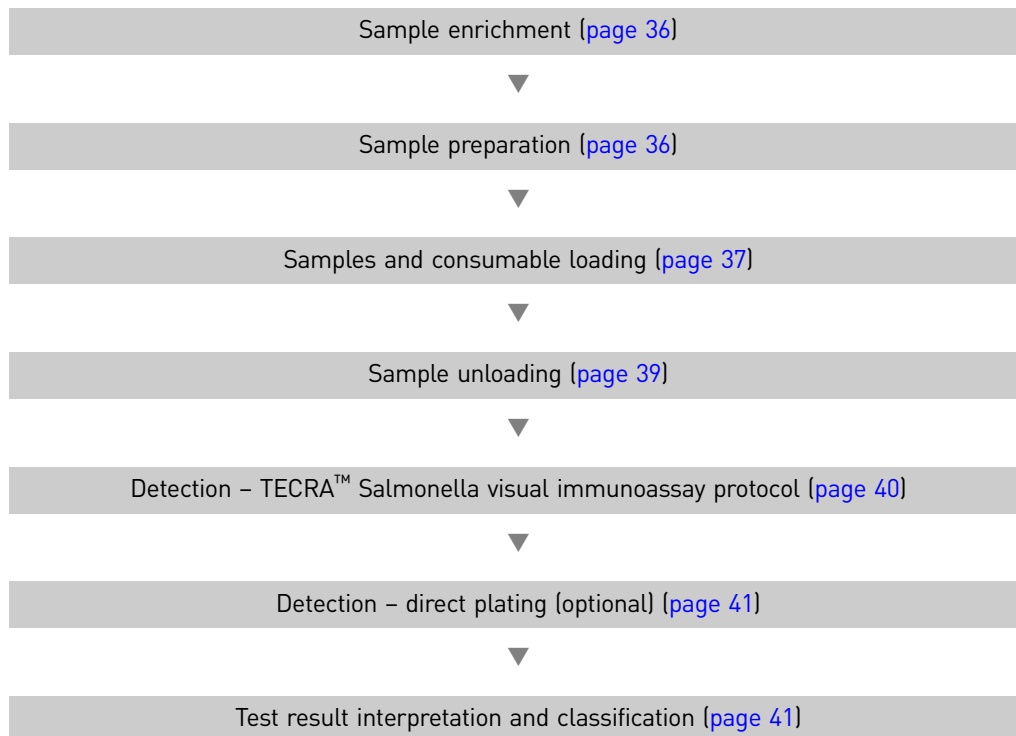

Procedural guidelines

- Use aseptic technique and good laboratory practices at all times.
- A facemask should be worn when weighing out powders.
- Care must be taken when boiling agar prior to autoclaving, and heat-resistant gloves should be worn when handling hot flasks of liquid.
- Take care when handling plates or vessels that contain microorganisms.
- Avoid generating aerosols, as pathogenic organisms may be present.
- Used or unused reagents, used media, sample enrichments, as well as other contaminated disposable materials should be disposed of following procedures for infectious or potentially infectious products.
- Sample enrichments might be contaminated with pathogenic organisms infectious to humans, so all waste must be treated as biohazardous and handled and disposed using safe laboratory practices, in accordance and compliance with all appropriate regulations.

Sample enrichment

Note: Certain food types and swabs/sponges can benefit from an alternative enrichment strategy (see [Appendix A, “Alternative Enrichment Methods”](#)).

1. Prepare a 1:10 dilution of the food sample in the appropriate **prewarmed** ($37 \pm 1^\circ\text{C}$) enrichment media in a sterile bag.

Note: For example, add 25 g of food sample to 225 mL of prewarmed media or add 325 g of food sample to 2925 mL of prewarmed media.

IMPORTANT! It is critical that the enrichment media is **prewarmed** to $37 \pm 1^\circ\text{C}$ prior to adding the food sample. To prevent cooling, all samples should then immediately be placed in a 37°C incubator.

2. Homogenize the sample by massaging the sample between the fingers for 10–15 seconds to disperse any large clumps of material (Narang and Cray, 2006).
3. Incubate at $37 \pm 1^\circ\text{C}$ for a **minimum** of **18–24 hours**.

Sample preparation

1. Remove the Sample Vessel and Elution Vessel from the consumable kit packaging and place into the Sample Vessel Holder.
2. Loosen the lids from the vessels and partially remove, leaving an opening through which to add sample and wash buffer.
3. Place 50 mL of your sample in the Sample Vessel.


Note: If the samples are highly particulate or has a high fat content, a Seward plain sterile bag with internal strainer may be used (Seward Product Code BA6041/STR).

4. Store the individual enriched samples at $5 \pm 3^\circ\text{C}$ for potential reanalysis until the test result is confirmed.


Note: Do not store for more than 32 hours. If refrigerated after pre-enrichment, samples should be rewarmed to $37 \pm 1^\circ\text{C}$ prior to removal of aliquots for analysis.

Samples and consumable loading


1. Add PBS (Cat. no. AM9624, diluted to 1X) to the fill line of the Elution Vessel (approximately 35 mL of 1X PBS).
2. Place the lids on the Sample and Elution Vessels, making sure that the vessels are sealed all the way around.
3. Ensure the Pathatrix® paramagnetic beads (included in the kit) are fully resuspended by agitating the bead vial (for example, by vortexing bead vial or inversion of sealed vial), and add 50 µL of the Pathatrix® paramagnetic bead suspension into the spout on the lid of the Sample Vessel.

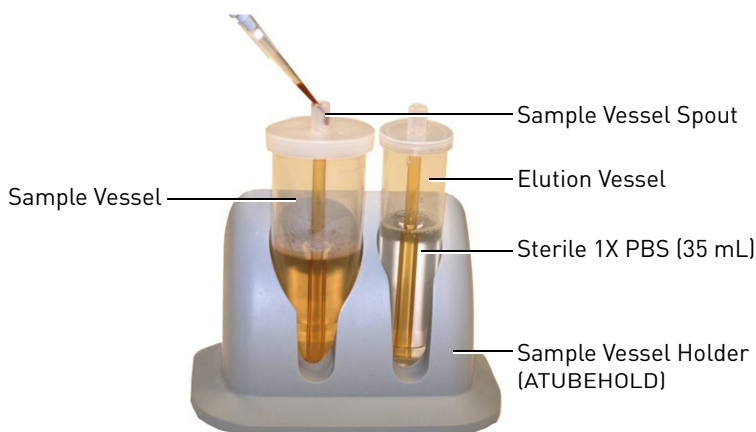

4. Remove the capture-phase kit from the bag and orient with the valve plunger pointing left. Connect the valve firmly to the lids of the Sample and Elution Vessels.

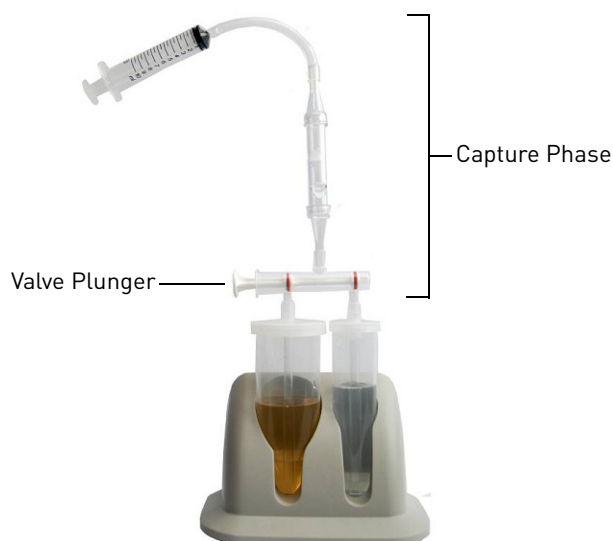

5. Holding both vessels, lift the assembled vessels and attached capture-phase kit out of the Sample Vessel Holder.
6. Place the vessels into the Cartridge, pushing them in firmly from the bottom upwards.
7. Firmly push the remainder of the kit into the Cartridge, starting with the valve followed by the capture phase and finally the syringe.
8. Push the magnet slider across into the locking position and press the magnet release button on the side of the Cartridge to ensure the kit is positioned correctly and the magnets will move away from the capture phase freely.

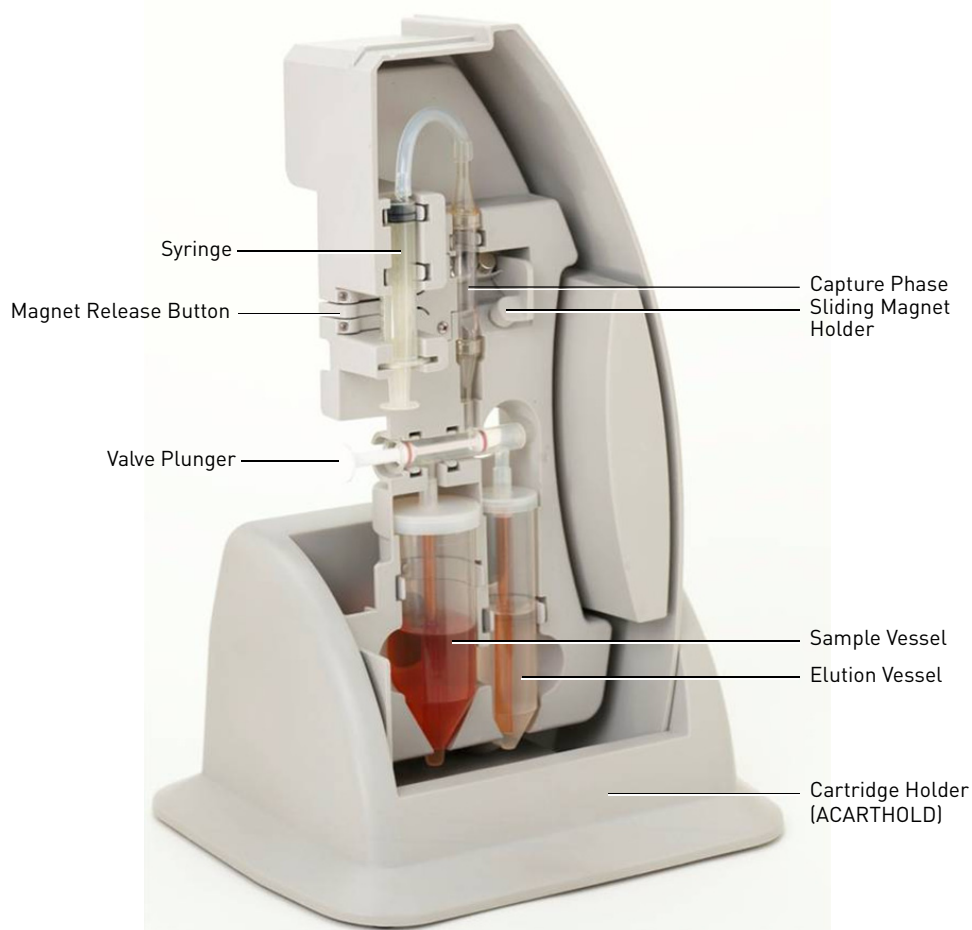

9. Reset the magnets into the locking position.
10. Insert the Cartridge into the Pathatrix® Auto Instrument until in the locking position.
11. Press the numbered button above the appropriate Cartridge. The LED above the Cartridge will turn green, indicating the run has begun (approximately 14 minutes).


Sample unloading

1. At the end of the run, the LED will flash red and green alternately.
2. Press the button above the appropriate Cartridge to initialize the draining step (approximately 1 minute).
3. When the draining step is complete and the LED is illuminated red, remove the Cartridge by pulling it out, away from the instrument.
4. Remove the syringe from the Cartridge and carefully pull out the rest of the kit, starting at the top and working downwards. When removing the Sample and Elution Vessels from the Cartridge, hold firmly by the vessels themselves to prevent spillage.
5. Place both vessels in the Tube Rack (also known as the Sample Vessel Holder or rack).
6. Remove the lid from the Elution Vessel and, leaving the Elution Vessel in place in the rack, lift away the rest of the consumable, including the Sample Vessel, and discard.
7. Place the Elution Vessel in the rack and cap with a flat lid (provided in the kit). Leave in the rack for 1 minute to allow capture of the Pathatrix® paramagnetic beads.
8. Remove all the liquid from the Elution Vessel, without removing the vessel from the rack, taking care not to disturb the captured Pathatrix® paramagnetic beads.
9. Remove the Elution Vessel from the vessel holder and use 500 µL of **prewarmed** BPW to resuspend the Pathatrix® paramagnetic beads.
10. Transfer 400 µL to a microcentrifuge tube.
Note: 100 µL of the Pathatrix® paramagnetic beads can be retained for direct plating prior to the additional 2-hour growth step.
11. Vortex the tube for 10 seconds.
12. Incubate the tube for 2 hours at 37 ±1°C.
13. Remove the tube from the incubator and vortex for 10 seconds.

Detection – TECRA™ *Salmonella* visual immunoassay protocol

IMPORTANT! It is solely the operator's responsibility to refer to the package insert instructions and instrument user manual for information regarding the correct set up and operation of the TECRA™ *Salmonella* visual immunoassay system.

1. Ensure that the tubes are securely capped and then heat in a water bath at 95°C for 15 minutes to release target antigens.
2. After heating, immediately place the tube into the DynaMag™-2 Magnet (Cat. no. 123.21D) for 2 minutes to allow the Pathatrix® paramagnetic beads to be pulled out of the lysate.
3. Once the Pathatrix® paramagnetic beads have aggregated against the magnet, open the tube and transfer 200 µL of the "bead-free" supernatant into the designated well of the TECRA™ microplate.
4. To proceed, refer to the TECRA™ operating instructions.

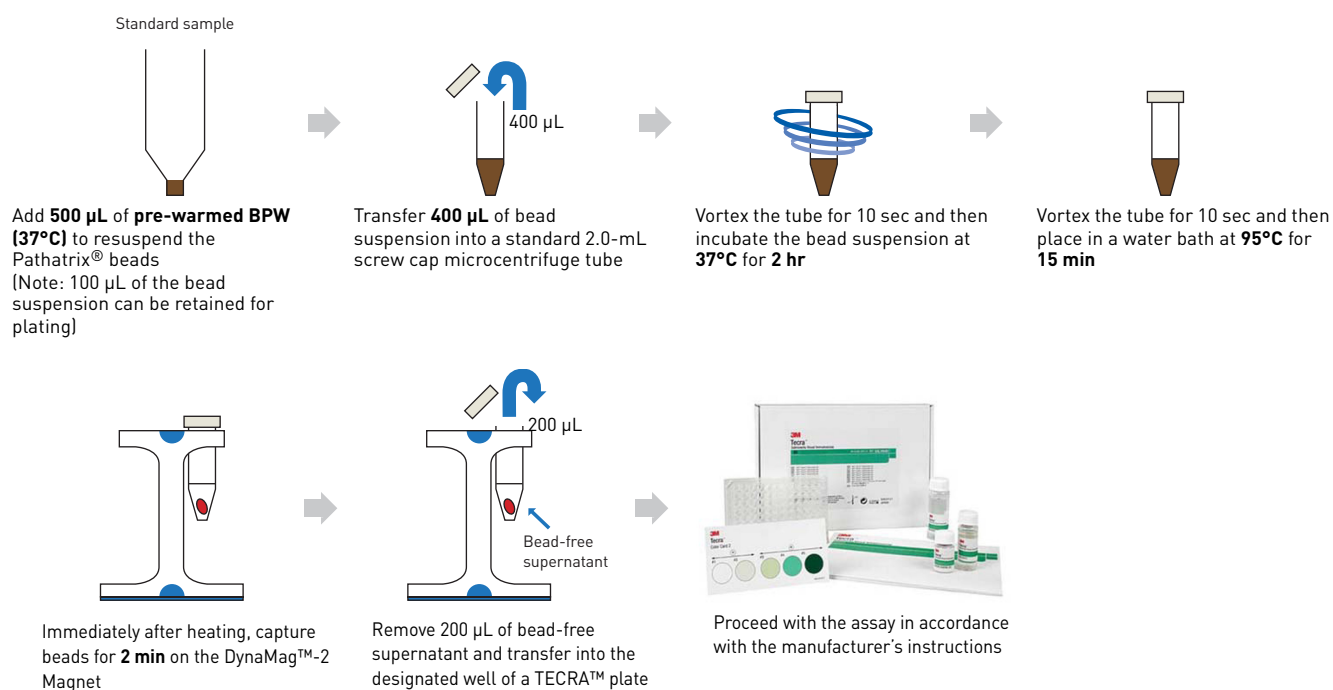

If a positive ELISA result is obtained, an aliquot of the Pathatrix® paramagnetic beads should be plated out (see the following section, "Detection – direct plating (optional)").

Detection – direct plating (optional)

Note: We recommend streaking the sample out to generate individual colonies, as opposed to spread plating.

1. Streak 50 µL of the unlysed Pathatrix® paramagnetic bead suspension over a well-dried XLD (xylose lysine desoxycholate) agar plate and another 50 µL onto an appropriate second selective plate medium.

Note: The laboratory may choose which medium to use, but the second selective plate should be any other solid selective medium complimentary to XLD and especially appropriate for the isolation of lactose-positive *Salmonella*, *Salmonella typhi*, and *Salmonella paratyphi* strains.

2. Allow the plates to dry for approximately 10 minutes then invert and incubate at the required temperature for 18–24 hours or as recommended by the manufacturer.

A presumptive positive result is defined as the isolation of typical, suspicious, or atypical *Salmonella* spp. colonies on the agar plates used. A presumptive positive isolate should be subsequently confirmed by the use of subculture and appropriate biochemical and serological tests as detailed in USDA Microbiology Laboratory Guidebook (MLG) 4.04 for cooked ham or FDA Bacteriological Analytical Manual (BAM), Chapter 5, sections D and E for tomatoes and chocolate (see [“References” on page 58](#)).

Test result interpretation and classification

The Pathatrix® *Salmonella* spp. Kit is designed as a sample preparation method for presence/absence detection of *Salmonella* spp. in food matrices

Using the Pathatrix® *Salmonella* spp. Kit linked to the TECRA™ *Salmonella* visual immunoassay, presumptive results can be obtained, prior to confirmation, within 22–28 hours.


Once confirmed, the results are reported as:

- *Salmonella* spp. **Detected** in 25–325 g (sample matrices)
- *Salmonella* spp. **Not detected** in 25–325 g (sample matrices)

5

Pathatrix® *Salmonella* spp. Kit Linked to the bioMérieux VIDAS® Immunoassay Testing System

Workflow

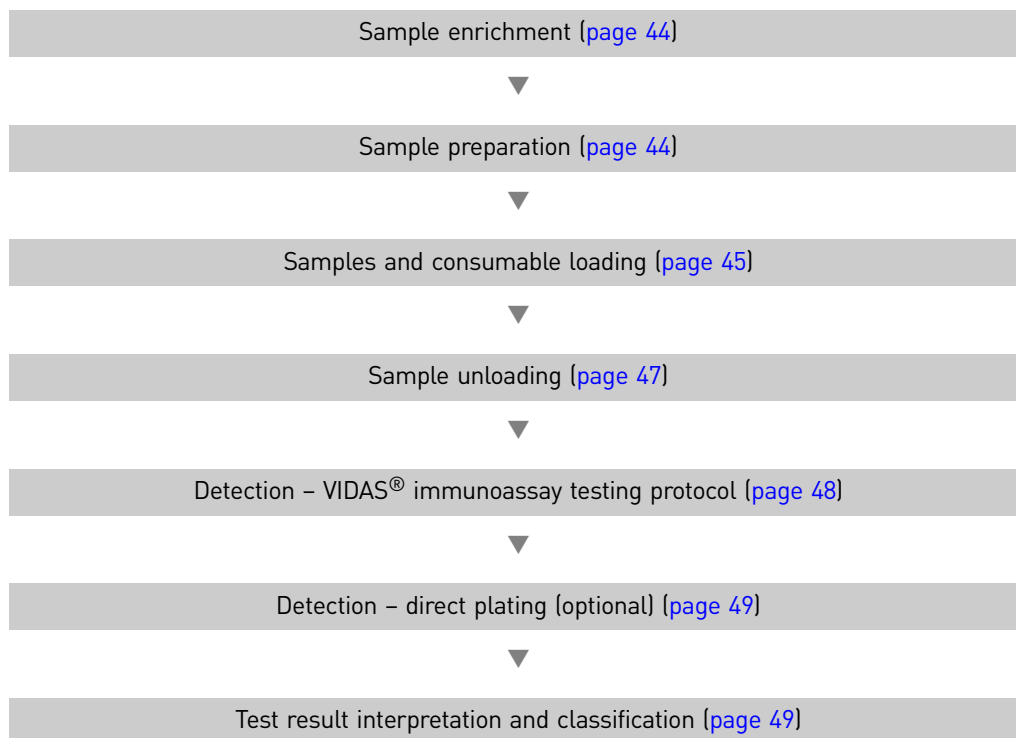

Procedural guidelines

- Use aseptic technique and good laboratory practices at all times.
- A facemask should be worn when weighing out powders.
- Care must be taken when boiling agar prior to autoclaving, and heat-resistant gloves should be worn when handling hot flasks of liquid.
- Take care when handling plates or vessels that contain microorganisms.
- Avoid generating aerosols, as pathogenic organisms may be present.
- Used or unused reagents, used media, sample enrichments, as well as other contaminated disposable materials should be disposed of following procedures for infectious or potentially infectious products.
- Sample enrichments might be contaminated with pathogenic organisms infectious to humans, so all waste must be treated as biohazardous and handled and disposed using safe laboratory practices, in accordance and compliance with all appropriate regulations.

Sample enrichment

Note: Certain food types and swabs/sponges can benefit from an alternative enrichment strategy (see [Appendix A, “Alternative Enrichment Methods”](#)).

1. Prepare a 1:10 dilution of the food sample in the appropriate **prewarmed** ($37 \pm 1^\circ\text{C}$) enrichment media in a sterile bag.

Note: For example, add 25 g of food sample to 225 mL of prewarmed media or add 325 g of food sample to 2925 mL of prewarmed media.

IMPORTANT! It is critical that the enrichment media is **prewarmed** to $37 \pm 1^\circ\text{C}$ prior to adding the food sample. To prevent cooling, all samples should then immediately be placed in a 37°C incubator.

2. Homogenize the sample by massaging the sample between the fingers for 10–15 seconds to disperse any large clumps of material (Narang and Cray, 2006).
3. Incubate at $37 \pm 1^\circ\text{C}$ for a **minimum** of **18–24 hours**.

Note: We recommend that these sub-samples for pooling are processed immediately or, if storage is required, that the samples are refrigerated at $5 \pm 3^\circ\text{C}$. Samples should be rewarmed to $37 \pm 1^\circ\text{C}$ prior to analysis on the Pathatrix® Auto Instrument. The remaining enriched sample should be stored until the results of the pooled sample have been determined.

Sample preparation

1. Remove the Sample Vessel and Elution Vessel from the consumable kit packaging and place into the Sample Vessel Holder.
2. Loosen the lids from the vessels and partially remove, leaving an opening through which to add sample and wash buffer.
3. Place 50 mL of your sample in the sample vessel.


Note: If the samples are highly particulate and/or contain a high fat content, a Seward plain sterile bag with internal strainer may be used (Seward Product Code BA6041/STR).

4. Store the individual enriched samples at $5 \pm 3^\circ\text{C}$ for potential reanalysis until the test result is confirmed.


Note: Do not store for more than 32 hours. If refrigerated at $5 \pm 3^\circ\text{C}$ prior to testing, samples should be rewarmed to $37 \pm 1^\circ\text{C}$ prior to removal of aliquots for analysis.

Samples and consumable loading


1. Add PBS (Cat. no. AM9624, diluted to 1X) to the fill line of the Elution Vessel (approximately 35 mL of 1X PBS).
2. Place the lids on the Sample and Elution Vessels, making sure that the vessels are sealed all the way around.
3. Ensure the Pathatrix® paramagnetic beads (included in the kit) are fully resuspended by agitating the bead vial (for example, by vortexing bead vial or inversion of sealed vial), and add 50 µL of the Pathatrix® paramagnetic bead suspension into the spout on the lid of the Sample Vessel.

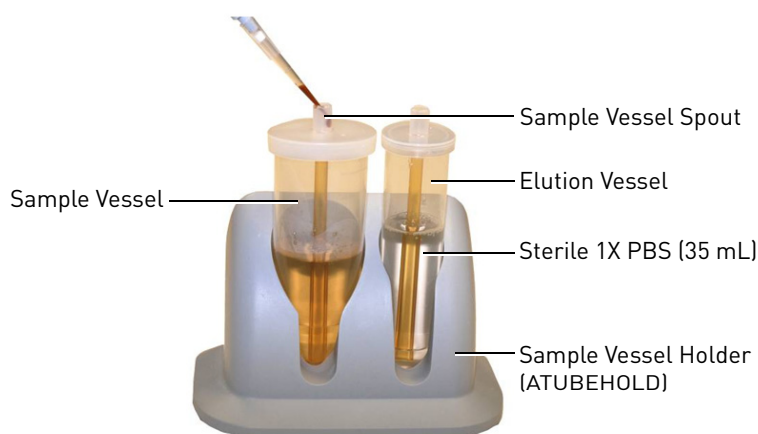

4. Remove the capture-phase kit from the bag and orient with the valve plunger pointing left. Connect the valve firmly to the lids of the Sample and Elution Vessels.

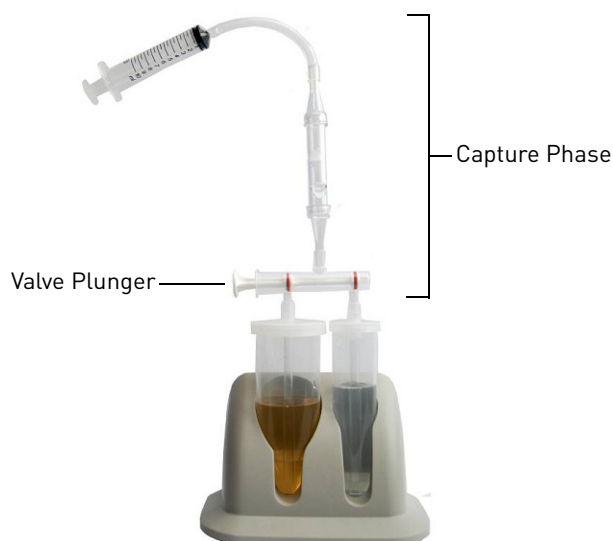

5. Holding both vessels, lift the assembled vessels and attached capture-phase kit out of the Sample Vessel Holder.
6. Place the vessels into the Cartridge, pushing them in firmly from the bottom upwards.
7. Firmly push the remainder of the kit into the Cartridge, starting with the valve, followed by the capture phase and finally by the syringe.
8. Push the magnet slider across into the locking position and press the magnet release button on the side of the Cartridge to ensure the kit is positioned correctly and the magnets will move away from the capture phase freely.

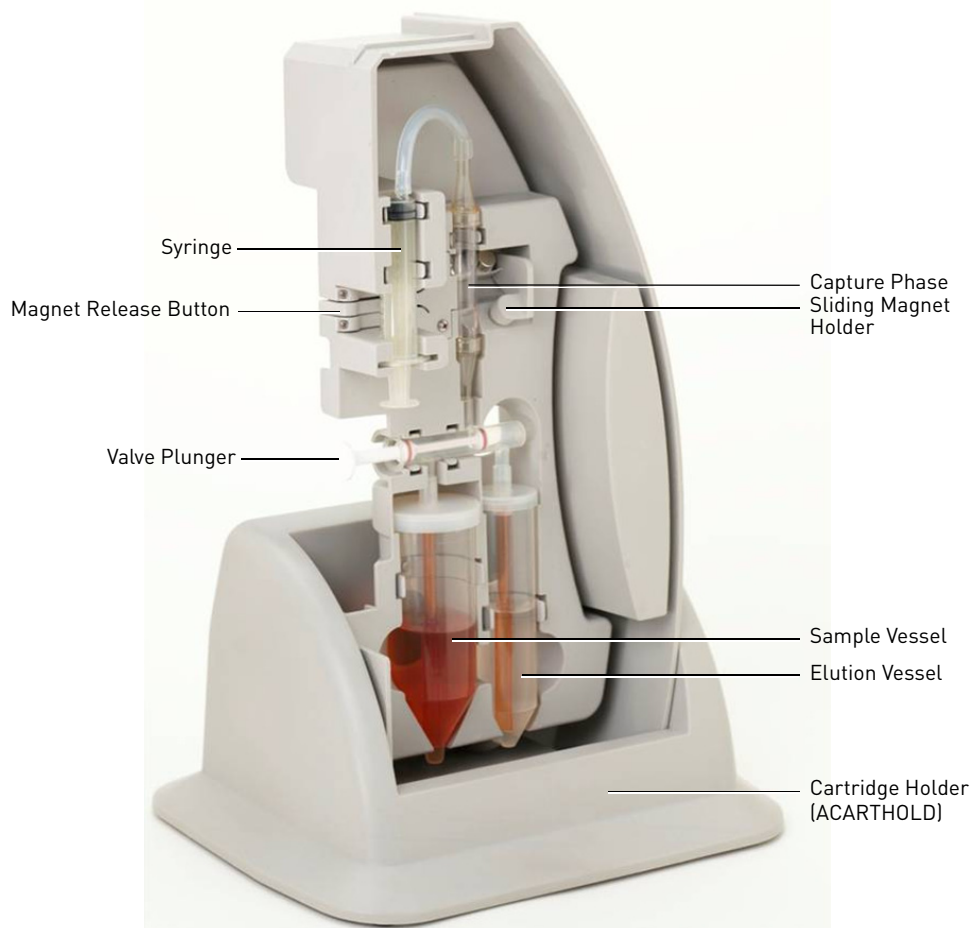

9. Reset the magnets into the locking position.
10. Insert the Cartridge into the Pathatrix® Auto Instrument until in the locking position.
11. Press the numbered button above the appropriate Cartridge. The LED above the Cartridge will turn green, indicating the run has begun (approximately 14 minutes).


Sample unloading

1. At the end of the run, the LED will flash red and green alternately.
2. Press the button above the appropriate Cartridge to initialize the draining step (approximately 1 minute).
3. When the draining step is complete and the LED is illuminated red, remove the Cartridge by pulling it out, away from the instrument.
4. Remove the syringe from the Cartridge and carefully pull out the rest of the kit, starting at the top and working downwards. When removing the Sample and Elution Vessels from the Cartridge, hold firmly by the vessels themselves to prevent spillage.
5. Place both vessels in the Tube Rack (also known as the Sample Vessel Holder or rack).
6. Remove the lid from the Elution Vessel and, leaving the Elution Vessel in place in the rack, lift away the rest of the consumable, including the Sample Vessel, and discard.
7. Place the Elution Vessel in the rack and cap with a flat lid (provided in the kit). Leave in the rack for 1 minute to allow capture of the Pathatrix® paramagnetic beads.
8. Remove all the liquid from the Elution Vessel, without removing the vessel from the rack, taking care not to disturb the captured Pathatrix® paramagnetic beads.
9. Remove the Elution Vessel from the vessel holder and use 700 µL of **prewarmed** BPW to resuspend the Pathatrix® paramagnetic beads.
10. Transfer 600 µL to a microcentrifuge tube.
Note: 100 µL of the Pathatrix® paramagnetic beads can be retained for direct plating prior to the additional 2-hour growth step.
11. Vortex the tube for 10 seconds.
12. Incubate the tube for 2 hours at 37 ±1°C.
13. Remove the tube from the incubator and vortex for 10 seconds.

Detection – VIDAS® immunoassay testing protocol

IMPORTANT! It is solely the operator's responsibility to refer to the package insert instructions and instrument user manual for information regarding the correct set up and operation of the VIDAS® system.

1. Ensure that the tubes are securely capped and then heat in a water bath at **95°C** for **15 minutes** to release target antigens.
2. After heating, immediately place the tube into the DynaMag™-2 Magnet (Cat. no. 123.21DK) for 1 minute to allow the Pathatrix® paramagnetic beads to be pulled out of the lysate.
3. Once the Pathatrix® paramagnetic beads have aggregated against the magnet, open the tube and remove 500 µL of the "bead-free" supernatant and inoculate the VIDAS® strip.
4. To proceed, please refer to the VIDAS® operating instructions.

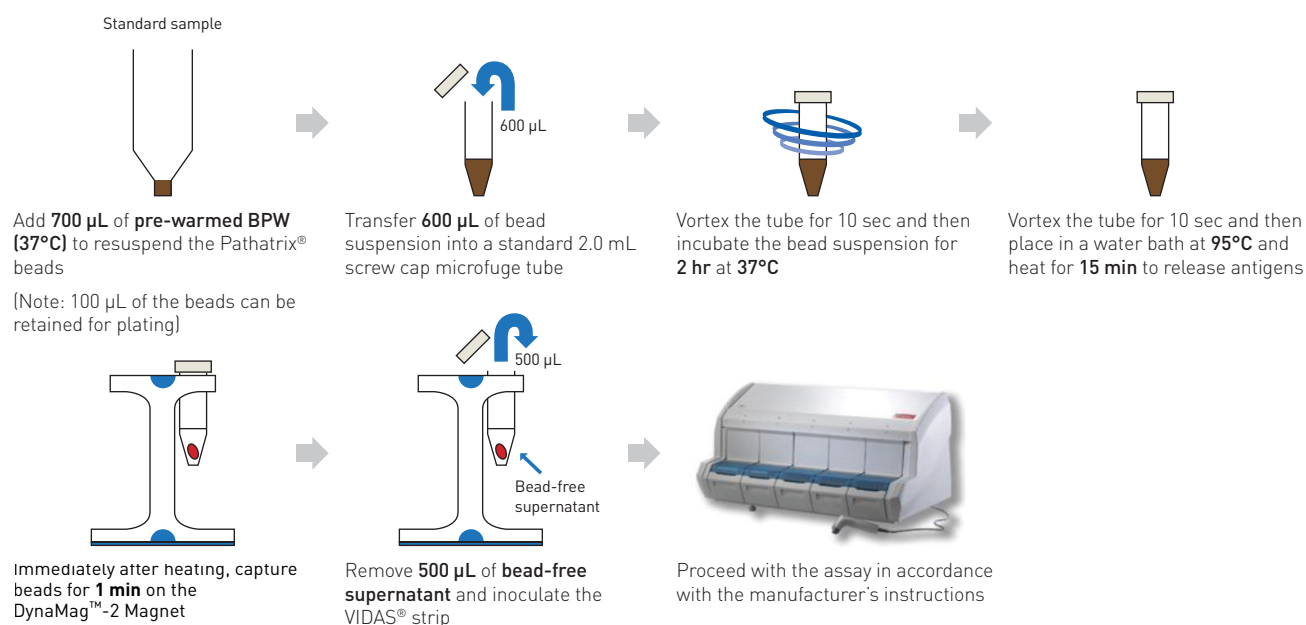

If a positive ELISA result is obtained, an aliquot of the Pathatrix® paramagnetic beads should be plated out (see the following section, "Detection – direct plating (optional)").

Detection – direct plating (optional)

Note: We recommend streaking the sample out to generate individual colonies, as opposed to spread plating.

1. Streak 50 µL of the unlysed Pathatrix® paramagnetic bead suspension over a well-dried XLD (xylose lysine desoxycholate) agar plate and another 50 µL onto an appropriate second selective plate medium.

Note: The laboratory may choose which medium to use, but the second selective plate should be any other solid selective medium complimentary to XLD and especially appropriate for the isolation of lactose-positive *Salmonella*, *Salmonella typhi*, and *Salmonella paratyphi* strains.

2. Allow the plates to dry for approximately 10 minutes then invert and incubate at the required temperature for 18–24 hours or as recommended by the manufacturer.

A presumptive positive result is defined as the isolation of typical, suspicious, or atypical *Salmonella* spp. colonies on the agar plates used. A presumptive positive isolate should be subsequently confirmed by the use of subculture and appropriate biochemical and serological tests as detailed in USDA Microbiology Laboratory Guidebook (MLG) 4.04 for cooked ham or FDA Bacteriological Analytical Manual (BAM), Chapter 5, sections D and E for tomatoes and chocolate (see [“References” on page 58](#)).

Test result interpretation and classification

The Pathatrix® *Salmonella* spp. Kit is designed as a sample preparation method for presence/absence detection of *Salmonella* spp. in food matrices

Using the Pathatrix® *Salmonella* spp. Kit linked to the VIDAS® immunoassay, presumptive results can be obtained, prior to confirmation, within 22–28 hours.

Once confirmed, the results are reported as:

- *Salmonella* spp. **Detected** in 25–325 g (sample matrices)
- *Salmonella* spp. **Not detected** in 25–325 g (sample matrices)

Alternative Enrichment Methods

Specific recommendations for high background microflora

For certain food groups, particularly where high levels of background microflora are present, the *Salmonella* assay results can be enhanced by the use of Tetrathionate (TT) broth instead of Buffered Peptone Water. The level of background contamination on the agar plates can be reduced by this method.

Examples of food groups that benefit from the use of TT broth are:

- Raw meat
- Raw vegetables
- Salads
- Fruits
- Ready meals

TT broth information

Materials

- Tetrathionate broth base is available from Oxoid (Product Code CM0029).
- Iodine Solution (to be added following preparation and heating of TT broth base).

Note: Prepare the iodine solution in advance to allow the iodine to dissolve, but use on the same day that it is made.

Directions for use

1. Add 77 g of TT broth base to 1 L of distilled water and bring to a boil.
2. Cool below 45°C.
3. Add 20 mL of iodine solution and mix well.

Following boiling, the prepared base can be stored for several weeks at $5 \pm 3^\circ\text{C}$, but once the iodine has been added, the media should be used the same day and any excess should be disposed of.

Specific recommendations for milk powder, chocolate, and cocoa-based samples

All milk powder, chocolate, or cocoa-based samples can benefit from an alternative enrichment.

Use prewarmed ($37 \pm 1^\circ\text{C}$) sterile UHT skim milk supplemented with Brilliant Green (0.002%), instead of Buffered Peptone Water, as the enrichment media.

- Milk powder samples should be incubated at $37 \pm 1^\circ\text{C}$ for 22 ± 2 hours.
- Chocolate and cocoa-based samples should be incubated, as described in the protocol, at $37 \pm 1^\circ\text{C}$ for 19 ± 1 hour.

Specific recommendations for potentially acidic/alkaline samples

Samples which potentially deviate from neutral pH should be prepared as follows:

1. Dilute the sample according to the sample enrichment protocol.
2. Incubate for 60 ± 5 minutes at room temperature.
3. Mix by hand massaging, and determine the pH.
4. If necessary, adjust the pH to 6.8 ± 0.2 , and mix well before determining the final pH.

Product overview

Description of target microorganisms

More than 2,400 *Salmonella* serotypes have been reported, all of which are potentially pathogenic. *Salmonella* is a frequently reported cause of foodborne illness, occurring in both epidemics and in isolated cases. *Salmonella* bacteria are the causative agent for Salmonellosis. Outbreaks have been associated with raw meats and poultry, eggs, milk and dairy products, seafood, coconut sauces, salad dressings, cocoa, chocolate, spices, frozen products, and vegetables such as hot peppers.

Audience

The Pathatrix® *Salmonella* spp. Kit is for professional use only and is intended for use by qualified users interested in determining the presence/absence of *Salmonella* spp. in food samples. Users may include, but are not limited to, food producers, food processors, food manufacturers, retailers, and microbiology testing laboratories.

Sampling protocol

The standard food sample size used in the Pathatrix® Auto system is 25 g of food diluted with 225 mL of enrichment medium. We recommend that the sub-samples removed for analysis are processed immediately or, if storage is required, that the samples are refrigerated at $5 \pm 3^{\circ}\text{C}$. Samples should be rewarmed to $37 \pm 1^{\circ}\text{C}$ prior to analysis with the Pathatrix® Auto system.

Kit sensitivity

The sample preparation procedure allows you to detect as few as 1–10 cfu from 25–325 g of food samples after enrichment. The limitation of the Pathatrix® *Salmonella* spp. Kit is in the ability of the target to reproduce in the enrichment medium, be captured by the magnet, and subsequently be detected by BAX® PCR, R.A.P.I.D® LT PCR, *Salmonella* OPTIMA, TECRA™ *Salmonella* visual immunoassay, or the VIDAS® immunoassay testing systems or be isolated on selective agar plates.

CAUTION! The Pathatrix® kit has been evaluated on cooked ham, chocolate, and chopped tomato food matrices. Given the wide variety of products and manufacturing procedures, we recommend that you check that the composition of the matrices to be tested does not affect the reliability of the results.

A negative result does not guarantee the absence of target organism in the original sample and may be due to the inability of the organism to adequately reproduce to required levels in the enrichment medium (with subsequent outgrowth on selective agar plates) potentially due to, but not limited to, competitive microflora, sub-lethal injury, or matrix inhibition.

Operating conditions

The Pathatrix® Auto Instrument is for indoor use only and for altitudes not exceeding 2000 m (6500 ft) above sea level.

Temperature and humidity requirements	
Condition	Acceptable range
Temperature	5–40°C
Humidity	Maximum relative humidity 80% for temperatures up to 31°C, decreasing to 50%

Related materials from Life Technologies

Item	Cat. no.
Related consumable kits with associated beads	
Pathatrix® <i>Salmonella</i> spp. Kit – Same Day	APS50SD
Pathatrix® 5-Pooling <i>Salmonella</i> spp. Kit	APS250P
Pathatrix® 10-Pooling <i>Salmonella</i> spp. Kit	APS500P
Pathatrix® 5-Pool DUAL (<i>E. coli</i> / <i>Salmonella</i> spp.) Kit	APDES250P
Pathatrix® DUAL (<i>Listeria</i> / <i>Salmonella</i> spp.) Kit	APD50
Equipment	
Pathatrix® Auto Instrument	PATHATRIXAUTO
Cartridge Rack (optional for use with the Pathatrix® Auto Instrument; holds 5 Cartridges)	ACARTRACK
DynaMag™-2 Magnet (for use with microcentrifuge tubes)	123.21D
Magnetic Capture Plate	MAGNETICPLATE
Reagents	
PBS, 10X, pH 7.4	AM9624 or AM9625
Related PCR assay	
MicroSEQ® <i>Salmonella</i> spp. Detection Kit	4403930



WARNING! GENERAL SAFETY. Using this product in a manner not specified in the user documentation may result in personal injury or damage to the instrument or device. Ensure that anyone using this product has received instructions in general safety practices for laboratories and the safety information provided in this document.

- Before handling chemicals, read and understand all applicable Safety Data Sheets (SDSs) and use appropriate personal protective equipment (gloves, gowns, eye protection, etc). To obtain SDSs, see the [“Documentation and Support”](#) section in this document.
-

Chemical safety



WARNING! GENERAL CHEMICAL HANDLING. To minimize hazards, ensure laboratory personnel read and practice the general safety guidelines for chemical usage, storage, and waste provided below, and consult the relevant SDS for specific precautions and instructions:

- Read and understand the Safety Data Sheets (SDSs) provided by the chemical manufacturer before you store, handle, or work with any chemicals or hazardous materials. To obtain SDSs, see the [“Documentation and Support”](#) section in this document.
- Minimize contact with chemicals. Wear appropriate personal protective equipment when handling chemicals (for example, safety glasses, gloves, or protective clothing).
- Minimize the inhalation of chemicals. Do not leave chemical containers open. Use only with adequate ventilation (for example, fume hood).
- Check regularly for chemical leaks or spills. If a leak or spill occurs, follow the manufacturer's cleanup procedures as recommended in the SDS.
- Handle chemical wastes in a fume hood.
- Ensure use of primary and secondary waste containers. (A primary waste container holds the immediate waste. A secondary container contains spills or leaks from the primary container. Both containers must be compatible with the waste material and meet federal, state, and local requirements for container storage.)
- After emptying a waste container, seal it with the cap provided.

- Characterize (by analysis if necessary) the waste generated by the particular applications, reagents, and substrates used in your laboratory.
- Ensure that the waste is stored, transferred, transported, and disposed of according to all local, state/provincial, and/or national regulations.
- **IMPORTANT!** Radioactive or biohazardous materials may require special handling, and disposal limitations may apply.

Specific chemical handling

CAS	Chemical	Phrase
26628-22-8	Sodium Azide	Sodium azide may react with lead and copper plumbing to form highly explosive metal azides.

Biological hazard safety



WARNING! BIOHAZARD. Biological samples such as tissues, body fluids, infectious agents, and blood of humans and other animals have the potential to transmit infectious diseases. Follow all applicable local, state/provincial, and/or national regulations. Wear appropriate protective equipment, which includes but is not limited to: protective eyewear, face shield, clothing/lab coat, and gloves. All work should be conducted in properly equipped facilities using the appropriate safety equipment (for example, physical containment devices). Individuals should be trained according to applicable regulatory and company/institution requirements before working with potentially infectious materials. Read and follow the applicable guidelines and/or regulatory requirements in the following:

In the U.S.:

- U.S. Department of Health and Human Services guidelines published in Biosafety in Microbiological and Biomedical Laboratories found at: www.cdc.gov/biosafety
- Occupational Safety and Health Standards, Bloodborne Pathogens (29 CFR§1910.1030), found at: www.access.gpo.gov/nara/cfr/waisidx_01/29cfr1910a_01.html
- Your company's/institution's Biosafety Program protocols for working with/handling potentially infectious materials.
- Additional information about biohazard guidelines is available at: www.cdc.gov

In the EU:

Check local guidelines and legislation on biohazard and biosafety precaution and refer to the best practices published in the World Health Organization (WHO) Laboratory Biosafety Manual, third edition, found at: www.who.int/csr/resources/publications/biosafety/WHO_CDS_CSR_LYO_2004_11/en/

Documentation and Support

Obtaining SDSs

Safety Data Sheets (SDSs) are available from www.lifetechnologies.com/support.

Note: For the SDSs of chemicals not distributed by Life Technologies, contact the chemical manufacturer.

Obtaining Certificates of Analysis

The Certificate of Analysis provides detailed quality control and product qualification information for each product. Certificates of Analysis are available on our website. Go to www.lifetechnologies.com/support and search for the Certificate of Analysis by product lot number, which is printed on the box.

Obtaining support

Support email: foodsafety@lifetech.com

For the latest services and support information for all locations, go to:

www.lifetechnologies.com/support

At the website, you can:

- Access worldwide telephone and fax numbers to contact Technical Support and Sales facilities
- Search through frequently asked questions (FAQs)
- Search for user documents, SDSs, application notes, formulations, handbooks, certificates of analysis, citations, and other product support documents
- Obtain information about customer training

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at www.lifetechnologies.com/termsandconditions. If you have any questions, please contact Life Technologies at www.lifetechnologies.com/support.

References

Food Safety and Inspection Service (USDA). 2008. Isolation And Identification of *Salmonella* From Meat, Poultry, Pasteurized Egg and Catfish Products. MLG 4.04. Microbiology Laboratory Guidebook.

Narang, N. and Cray, W.C. 2006. Evaluation of Hand Mixing of Ground Beef and Poultry Samples as an Alternative to Stomaching for the Detection of *Salmonella*. *Food Protection Trends*. 26:14–19.

US FDA Bacteriological Analytical Manual (BAM), Chapter 5; go to www.fda.gov/Food/ScienceResearch/LaboratoryMethods/BacteriologicalAnalyticalManualBAM/default.htm and scroll to *Salmonella*.

Headquarters

5791 Van Allen Way | Carlsbad, CA 92008 USA | Phone +1 760 603 7200 | Toll Free in USA 800 955 6288

For support visit www.lifetechnologies.com/support or email techsupport@lifetech.com

www.lifetechnologies.com

October 2012

